



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
1-13-2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/7/2016 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 12/10/2015 and 1/7/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


1/12/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 12th, day of January, 2016.

Elma Escarcega

Notary Public

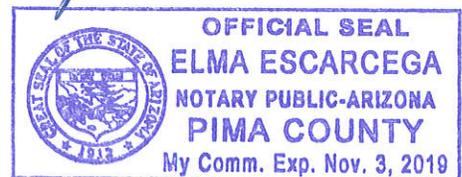

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TSK 11-13-2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 1/7/2016 TIME 15:36

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE William Loomis

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
01/07/2016 15:38:38

2 of 7
TSG WLL

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:39:30
Diagnostic Test	Pass	15:40:05
Air Blank	0.000	15:40:32
0.100 Cal Check	0.096	15:40:53
Air Blank	0.000	15:41:23
Subject Test	0.000	15:41:42
Air Blank	0.000	15:42:09
Wait	ABT*	
Air Blank	0.000	15:42:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
01/07/2016 15:42:48

3 of 7
WLL
T04

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:43:38
Diagnostic Test	Pass	15:44:13
Air Blank	0.000	15:44:41
0.100 Cal Check	RFI*	15:44:56
Air Blank	IPS**	15:45:03

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
01/07/2016 15:45:12

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:46:05
Diagnostic Test	Pass	15:46:40
Air Blank	0.000	15:47:07
0.100 Cal Check	0.095	15:47:28
Air Blank	0.000	15:47:57
Subject Test	0.110*	15:51:05
Air Blank	0.000	15:51:37
Air Blank	0.000	15:52:04
Subject Test	XXX**	15:52:25
Air Blank	0.000	15:52:59

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
WLL
TSG

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
01/07/2016 15:53:48

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:54:18
0.100 Cal Check	0.094	15:54:39
Air Blank	0.000	15:55:09

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLV
BL

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TS4
1.13.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 1/7/2016 TIME 15:36
Intoxilyzer Serial # 80-001163 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
01/07/2016 15:36:55

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:37:25
0.100 Cal Check	0.095	15:37:46
Air Blank	0.000	15:38:15

Preventative Maintenance
Performed -- Timer Reset

7 of 7
0 WLL
TSh