



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 4.7-2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/31/2016 Cal 1 0.095 AC Cal 2 0.095 AC

Between the dates of 3/3/2016 and 3/31/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
 The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


4/4/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 4th, day of April, 2016.

Elma Escarcega

Notary Public

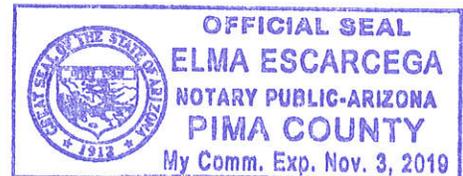

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
B4 4.7.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 3/31/2016 TIME 15:36

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE William Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
03/31/2016

15:38:35

Standard Lot#: 11514100A2  
Last Changed By: J. TRAINING #99999

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:40:14
Diagnostic Test	Pass	15:40:49
Air Blank	0.000	15:41:16
0.100 Cal Check	0.095	15:41:37
Air Blank	0.000	15:42:06
Subject Test	0.000	15:42:33
Air Blank	0.000	15:43:00
Wait	ABT*	
Air Blank	0.000	15:43:36

\*Sequence Aborted

2 of 7  
WLL  
TSH

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
03/31/2016 15:43:49

3 of 7  
WVW  
TSh

Standard Lot#: 11514100A2  
Last Changed By: J. TRAINING #99999

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:45:32
Diagnostic Test	Pass	15:46:07
Air Blank	0.000	15:46:34
0.100 Cal Check	RFI*	15:46:47
Air Blank	IPS**	15:46:58

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
03/31/2016

15:47:09

4 of 7  
WLV  
TSL

Standard Lot#: 11514100A2  
Last Changed By: J. TRAINING #99999

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:48:26
Diagnostic Test	Pass	15:49:01
Air Blank	0.000	15:49:29
0.100 Cal Check	0.095	15:49:49
Air Blank	0.000	15:50:19
Subject Test	0.349*	15:53:26
Air Blank	0.000	15:53:59
Air Blank	0.000	15:54:26
Subject Test	XXX**	15:54:58
Air Blank	0.000	15:55:31

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
03/31/2016 15:56:11

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. TRAINING #99999  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:56:40
0.100 Cal Check	0.095	15:57:01
Air Blank	0.000	15:57:30

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
6 WLL  
TJH

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TSJ  
4-7-2016  
WLL

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 3/31/2016 TIME 15:36  
Intoxilyzer Serial # 80-001163 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 095 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
03/31/2016 15:36:37

7 of 7  
TSL  
WLV

-----  
Standard Lot#: 11514100A2  
Last Changed By: J. TRAINING #99999  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:37:15
0.100 Cal Check	0.095	15:37:36
Air Blank	0.000	15:38:05

Preventative Maintenance  
Performed -- Timer Reset