



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9.22.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-16-16 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 8-18-16 and 9-16-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Amy L. Autret

Amy L. Autret


9-19-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 19th, day of September, 2016.

Elma Escarcega

Notary Public

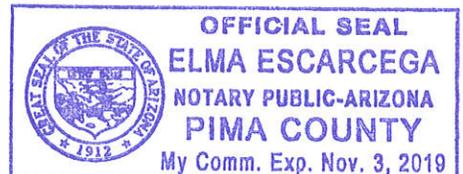

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TSK 9.22.2016

QA SPECIALIST Amy Rutz + 53546 AGENCY Tucson Police Department

DATE 9-16-16 TIME 1309

INTOXILYZER SERIAL # 80-001163 LOCATION Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

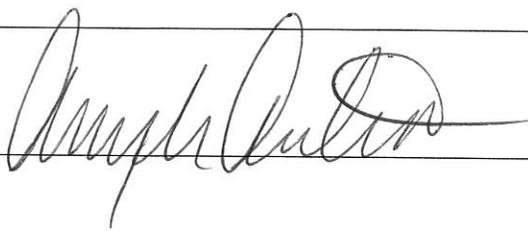SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
09/16/2016 13:11:27

2 of 7
W
Bh

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:12:56 |
| Diagnostic Test | Pass | 13:13:31 |
| Air Blank | 0.000 | 13:13:58 |
| 0.100 Cal Check | 0.096 | 13:14:19 |
| Air Blank | 0.000 | 13:14:49 |
| Subject Test | 0.000 | 13:15:17 |
| Air Blank | 0.000 | 13:15:45 |
| Wait | ABT* | |
| Air Blank | 0.000 | 13:16:22 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
09/16/2016 13:16:32

3 of 7
[Signature]
734

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:18:18 |
| Diagnostic Test | Pass | 13:18:53 |
| Air Blank | 0.000 | 13:19:20 |
| 0.100 Cal Check | RFI* | 13:19:33 |
| Air Blank | IPS** | 13:19:44 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
09/16/2016 13:20:02

4 of 7
DWB
134

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:21:04 |
| Diagnostic Test | Pass | 13:21:39 |
| Air Blank | 0.000 | 13:22:07 |
| 0.100 Cal Check | 0.096 | 13:22:27 |
| Air Blank | 0.000 | 13:22:57 |
| Subject Test | 0.293* | 13:26:04 |
| Air Blank | 0.000 | 13:26:37 |
| Air Blank | 0.000 | 13:27:05 |
| Subject Test | XXX** | 13:27:27 |
| Air Blank | 0.000 | 13:28:02 |

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
09/16/2016 13:29:17

5 of 7
JLD
TSh

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:29:46 |
| 0.100 Cal Check | 0.095 | 13:30:07 |
| Air Blank | 0.000 | 13:30:37 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000


STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TJH
9.22.2016

QA Specialist Amy Antret 53546 AGENCY Tucson Police Department
DATE 9-16-16 TIME 1309
Intoxilyzer Serial # 80-001163 Location Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

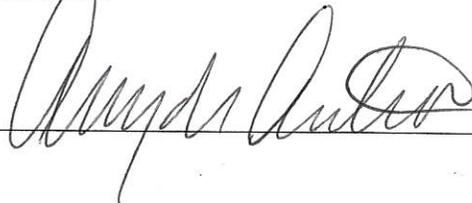

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
09/16/2016 13:09:46

7057
JL
TS4

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:10:18 |
| 0.100 Cal Check | 0.096 | 13:10:39 |
| Air Blank | 0.000 | 13:11:09 |

Preventative Maintenance
Performed -- Timer Reset