



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

je

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12.22.2006 Cal 1 0.097 AC Cal 2 0.098 AC

Between the dates of 12.1.2006 and 12.22.2006 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .


Terry S. Gallegos

Terry S. Gallegos

1.4.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

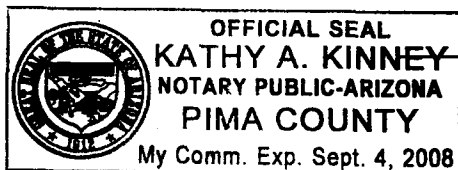

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 20942 AGENCY Tucson Police Department
DATE 12.22.2006 TIME 1656
INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9
DL

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *mu*
Serial Number : 80-001165
12/22/2006 16:57:58

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:58:52
Diagnostic Test	Pass	16:59:20
Air Blank	0.000	16:59:47
0.100 Cal Check	0.097	17:00:06
Air Blank	0.000	17:00:35
Subject Test	0.000	17:00:53
Air Blank	0.000	17:01:20
Five Minute Wait	ABT*	
Air Blank	0.000	17:01:55

*Sequence Aborted

Not a Successfully Completed Test Sequence

3079
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
12/22/2006

17:02:12

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:03:13
Diagnostic Test	Pass	17:03:40
Air Blank	IPS*	17:03:48
Air Blank	0.000	17:04:17

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
12/22/2006

17:19:06

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:20:00
Diagnostic Test	Fail*	17:20:20
Air Blank	RFI**	17:20:42

* Analytical Stability Test
**RFI Detect

Not a Successfully Completed Test Sequence

5079
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT *je*
Serial Number : 80-001165
12/22/2006 17:04:39

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:05:36
Diagnostic Test	Pass	17:06:03
Air Blank	0.000	17:06:31
0.100 Cal Check	0.098	17:06:50
Air Blank	0.000	17:07:19
Subject Test	0.201*	17:10:28
Air Blank	0.000	17:11:00
Air Blank	0.000	17:11:27
Subject Test	0.076	17:11:48
Air Blank	0.000	17:12:22
Five Minute Wait	ABT**	
Air Blank	0.000	17:13:04

*Deficient Sample

**Sequence Aborted

Invalid sample test accepted as 0.076 result

Not a Successfully Completed Test Sequence

6029
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
12/22/2006 17:13:26

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:14:20
Diagnostic Test	Pass	17:14:48
Air Blank	0.000	17:15:15
0.100 Cal Check	0.097	17:15:34
Air Blank	0.000	17:16:04
Subject Test	XXX*	17:16:56
Air Blank	0.000	17:17:30

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

12/22/2006

17:21:37

Standard Lot#: 502602I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:22:05
Cal Check	0.098	17:22:23
Air Blank	0.000	17:22:53


Operator's Signature

EXHIBIT G-2

8079
Ph

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

je

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28947 AGENCY Tucson Police Department
DATE 12-22-2006 TIME 1656
Intoxilyzer Serial # 80-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

12/22/2006

16:56:22


Standard Lot#: 502602I

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:56:50
Cal Check	0.097	16:57:08
Air Blank	0.000	16:57:38


Operator's Signature