



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

By
12-13-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-23-07 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 10-25-07 and 11-23-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

12-7-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 7, day of December, 2007.

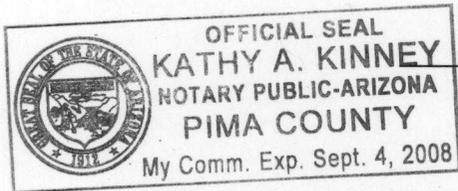

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1067 cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

154
12-13-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 11/23/07 TIME 08:15

INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

2 of 7 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
11/23/2007

08:19:32

By

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:20:28
Diagnostic Test	Pass	08:20:55
Air Blank	0.000	08:21:23
0.100 Cal Check	0.099	08:21:42
Air Blank	0.000	08:22:11
Subject Test	0.000	08:22:28
Air Blank	ABT*	08:22:38
Air Blank	0.000	08:23:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 26-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3067 cm

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
11/23/2007

08:23:30

BH

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:24:05
Air Blank	IPS**	08:24:13

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
11/23/2007 08:25:00

134

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:25:49
Diagnostic Test	Pass	08:26:17
Air Blank	0.000	08:26:44
0.100 Cal Check	0.100	08:27:03
Air Blank	0.000	08:27:33
Subject Test	0.133*	08:30:41
Air Blank	0.000	08:31:13
Air Blank	0.000	08:31:40
Subject Test	XXX**	08:32:02
Air Blank	0.000	08:32:36

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 26-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
11/23/2007 08:33:16

5-107 cu
B4

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:33:43
Cal Check	0.098	08:34:02
Air Blank	0.000	08:34:32

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

b507cu
154

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 11/23/07 TIME 08:15
Intoxilyzer Serial # S2-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

767cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
11/23/2007 08:17:24

B4

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:17:51
Cal Check	0.099	08:18:10
Air Blank	0.000	08:18:40

Clare Conley
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1062 cu

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4 12.13.07

Date: November 23, 2007

QAS: Clare Conley

Intoxilyzer Serial Number: 80-001165

New Lot Number: 725402I

New Expiration Date: September 12, 2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

2072 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
11/23/2007 08:15:59

TS4

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:16:27
Cal Check	0.099	08:16:45
Air Blank	0.000	08:17:15

Clare Conley
Operator's Signature