



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-001165 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

TS4
3.17.08

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-7-08 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2-8-08 and 3-7-08
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

3-14-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 14, day of
March, 2008.

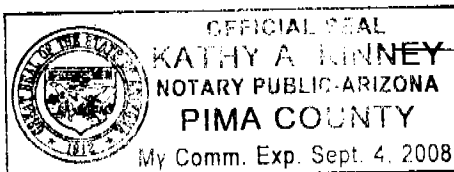

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the
information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

157
cu TB4
3.17.08QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 3-7-08 TIME 0739
INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 099
- AC

Instrument is operating properly and accurately. YES ☒ NO ☐COMMENTS NoneSIGNATURE Clare Conley

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
03/07/2008
07:42:17

Standard Lot#: 7254021
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M
Weight: 111
15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:43:08
Diagnostic Test	Pass	07:43:36
Air Blank	0.000	07:44:03
0.100 Cal Check	0.100	07:44:22
Air Blank	0.000	07:44:51
Subject Test	0.000	07:45:07
Air Blank	ABT*	07:45:14
Air Blank	0.000	07:45:42

*Sequence Aborted

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

03/07/2008

07:45:47

Standard Lot#: 725402I

Last Changed By: C. CONLEY

#34517

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:46:21
Air Blank	IPS**	07:46:29

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

03/07/2008

07:46:50

Standard Lot#: 725402I

Last Changed By: C. CONLEY

#34517

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:47:43
Diagnostic Test	Pass	07:48:11
Air Blank	0.000	07:48:38
0.100 Cal Check	0.100	07:48:57
Air Blank	0.000	07:49:27
Subject Test	0.067*	07:52:35
Air Blank	0.000	07:53:07
Air Blank	0.000	07:53:35
Subject Test	XXX**	07:53:57
Air Blank	0.000	07:54:30

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

03/07/2008

07:54:48


Standard Lot#: 725402I

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:55:15
Cal Check	0.099	07:55:34
Air Blank	0.000	07:56:04

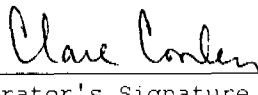

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3 7 08 TIME 0739
Intoxilyzer Serial # 50-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

03/07/2008

07:39:14

Standard Lot#: 725402I

Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:39:41
Cal Check	0.100	07:40:00
Air Blank	0.000	07:40:30

Clare Conley

Operator's Signature

787
Cm
134