



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
6-30-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-27-08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 5-30-08 and 6-27-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

6-27-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 27, day of June, 2008.

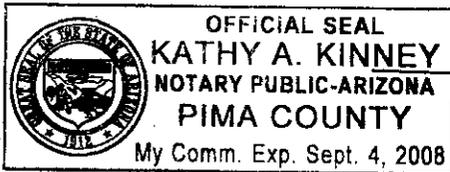

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
Ch
754
6-30-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 6-27-08 TIME 0733
INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(/) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Added six minutes to time on clock.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
06/27/2008 07:34:31

2 of 7
CW

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

104

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:35:26
Diagnostic Test	Pass	07:35:53
Air Blank	0.000	07:36:21
0.100 Cal Check	0.098	07:36:39
Air Blank	0.000	07:37:10
Subject Test	0.000	07:37:35
Air Blank	ABT*	07:37:42
Air Blank	0.000	07:38:10

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
06/27/2008 07:44:22

3 of 7
em

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

TS4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:45:04
Air Blank	IPS**	07:45:17

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
06/27/2008

07:45:23

4567
ew

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

134

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:46:13
Diagnostic Test	Pass	07:46:41
Air Blank	0.000	07:47:08
0.100 Cal Check	0.099	07:47:27
Air Blank	0.000	07:47:57
Subject Test	0.131*	07:51:06
Air Blank	0.000	07:51:39
Air Blank	0.000	07:52:06
Subject Test	XXX**	07:52:31
Air Blank	0.000	07:53:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
06/27/2008 07:53:46

5367
cu

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

TS4

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:54:13
Cal Check	0.098	07:54:32
Air Blank	0.000	07:55:02

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6/27
cu
TSY
6-30-08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-27-08 TIME 0733
Intoxilyzer Serial # 80-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
06/27/2008 07:33:08

7867
cu

Standard Lot#: 725402I
Last Changed By: C. CONLEY #34517

BY

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:33:35
Cal Check	0.099	07:33:54
Air Blank	0.000	07:34:24

Clare Conley

Operator's Signature