



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

*TSY*  
*5.11.2009*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/1/09 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 4/2/09 and 5/1/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
 The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


5/4/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 4, day of May, 2009.

A.M. Hernandez

Notary Public

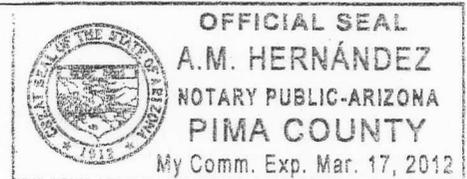

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 sc  
TSC  
5-11-2009

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 5/1/09 TIME 10:53  
INTOXILYZER SERIAL # 80-001165 LOCATION Team 1

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7 sc

124

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
05/01/2009 10:54:43

Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:55:25
Diagnostic Test	Pass	10:55:53
Air Blank	0.000	10:56:20
0.100 Cal Check	0.100	10:56:39
Air Blank	0.000	10:57:08
Subject Test	0.000	10:57:36
Air Blank	0.000	10:58:04
Five Minute Wait	ABT*	
Air Blank	0.000	10:58:38

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc  
184

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
05/01/2009 10:58:42

---

Standard Lot#: 804302I  
Last Changed By: S. CIONE - #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:59:13
Air Blank	IPS**	10:59:21

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
05/01/2009

10:59:26

4 of 7<sub>sc</sub>

154

Standard Lot#: 804302I  
Last Changed By: S. CIONE -#50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:00:11
Diagnostic Test	Pass	11:00:38
Air Blank	0.000	11:01:06
0.100 Cal Check	0.099	11:01:24
Air Blank	0.000	11:01:55
Subject Test	0.100*	11:05:03
Air Blank	0.000	11:05:34
Air Blank	0.000	11:06:02
Subject Test	XXX**	11:06:19
Air Blank	0.000	11:06:53

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
05/01/2009

11:07:12

5 of 7 sc

134

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:07:39
Cal Check	0.099	11:07:58
Air Blank	0.000	11:08:28

Sarah Cione  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

T34  
5-11-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/1/09 TIME 10:53  
Intoxilyzer Serial # 80-001165 Location Team 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
05/01/2009 10:53:10

7 of 7 sc

154

-----  
Standard Lot#: 804302I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:53:38
Cal Check	0.099	10:53:57
Air Blank	0.000	10:54:26

*Sarah Cione*

Operator's Signature