



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

754  
7-28-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/27/09 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 06/26/09 and 07/27/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

*Cheryl K. Straub*

Cheryl K. Straub


07/28/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 28, day of July, 2009.

*A.M. Hernandez*

Notary Public

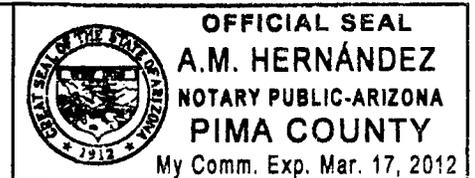

\*Includes documents specified in A.A.C. R13-10-104.A.5

199  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

154  
7-28-09

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/27/09 TIME 15:43

INTOXILYZER SERIAL # 80-001165 LOCATION TEAM I

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

Ambient Failure of First AIR blank  
After last <sup>First</sup> Calibration Check. Mouth spray  
sprayed in close proximity of BREATH HOSE.  
Calibration check was completed before this failure

SIGNATURE

Cheryl Straub

289  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
07/27/2009 15:45:24

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

BY

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:46:20
Diagnostic Test	Pass	15:46:47
Air Blank	0.000	15:47:14
0.100 Cal Check	0.098	15:47:33
Air Blank	0.000	15:48:03
Subject Test	0.000	15:48:25
Air Blank	ABT*	15:48:31
Air Blank	0.000	15:49:00

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

379  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
07/27/2009 15:55:44

TS4

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:56:20
Air Blank	IPS**	15:56:29

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

489  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
07/27/2009 15:49:21

TS4

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:50:13
Diagnostic Test	Pass	15:50:41
Air Blank	0.000	15:51:08
0.100 Cal Check	0.098	15:51:27
Air Blank	0.000	15:51:57
Subject Test	XXX*	15:54:52
Air Blank	0.000	15:55:25

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

JCS  
53?  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
07/27/2009 15:58:44

154

Standard Lot#: 826001I  
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:59:39
Diagnostic Test	Pass	16:00:06
Air Blank	0.000	16:00:34
0.100 Cal Check	0.098	16:00:53
Air Blank	0.000	16:01:23
Subject Test	0.282*	16:04:31
Air Blank	0.000	16:05:04
Air Blank	ABT**	16:05:13
Air Blank	0.000	16:05:41

\*Deficient Sample  
\*\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

679  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/27/2009

16:06:16

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

BY

-----  
QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:06:43
Cal Check	0.097	16:07:02
Air Blank	0.000	16:07:32

  
Operator's Signature

799  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

154  
7-28-09

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/27/09 TIME 15:43  
Intoxilyzer Serial # 80-001165 Location TEAM I

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

899  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/27/2009

15:43:24

AK

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

-----  
QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:43:52
Cal Check	0.097	15:44:10
Air Blank	0.000	15:44:41



Operator's Signature

999  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/27/2009

15:57:41

154

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	AMB*	15:58:09
Air Blank	0.000	15:58:37

\*Ambient Fail



\_\_\_\_\_  
Operator's Signature