



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
8.31.2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 08-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/27/09 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 07/27/09 and 08/27/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 20 09.

A.M. Hernández

Notary Public

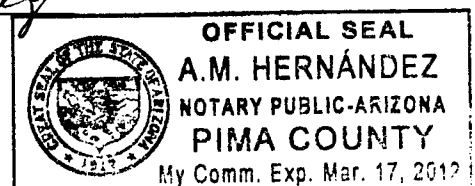

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

197
CS

B4
8.31.09

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department
DATE 08/27/09 TIME 07:10
INTOXILYZER SERIAL # 80-001165 LOCATION TEAM I

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

287
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
08/27/2009

07:12:45

154

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TPD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:13:53
Diagnostic Test	Pass	07:14:21
Air Blank	0.000	07:14:48
0.100 Cal Check	0.099	07:15:07
Air Blank	0.000	07:15:37
Subject Test	0.000	07:15:59
Air Blank	ABT*	07:16:05
Air Blank	0.000	07:16:33

*Sequence Aborted

Not a Successfully Completed Test Sequence

387
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
08/27/2009

07:17:01

134

Standard Lot#: 826001I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:17:51
Air Blank	IPS**	07:17:59

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

08/27/2009

07:18:20

by

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:19:15
Diagnostic Test	Pass	07:19:43
Air Blank	0.000	07:20:10
0.100 Cal Check	0.098	07:20:29
Air Blank	0.000	07:20:59
Subject Test	0.108*	07:24:08
Air Blank	0.000	07:24:40
Air Blank	0.000	07:25:07
Subject Test	XXX**	07:25:35
Air Blank	0.000	07:26:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

597
cs

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

08/27/2009

07:26:54

BY


Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:27:22
Cal Check	0.098	07:27:41
Air Blank	0.000	07:28:10

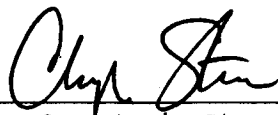

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 08/27/09 TIME 07:10
Intoxilyzer Serial # 80-001165 Location TEAM I

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

08/27/2009

07:10:35

Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

124

Test	g/210L	Time
Air Blank	0.000	07:11:03
Cal Check	0.099	07:11:21
Air Blank	0.000	07:11:52



Operator's Signature