



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/11/10 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 02/11/10 and 03/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

03/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th, day of March, 2010.

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

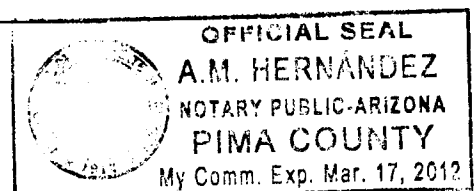


EXHIBIT G-4

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
3-11-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 03/11/10 TIME 07:55
INTOXILYZER SERIAL # 90-001165 LOCATION TEAM 1

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

277
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

03/11/2010

07:57:35

TS4

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:58:28
Diagnostic Test	Pass	07:58:56
Air Blank	0.000	07:59:23
0.100 Cal Check	0.098	07:59:42
Air Blank	0.000	08:00:12
Subject Test	0.000	08:00:30
Air Blank	ABT*	08:00:37
Air Blank	0.000	08:01:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

03/11/2010

08:01:08

TS4

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:01:42
Air Blank	IPS**	08:01:50

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

03/11/2010

08:01:59

154

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:02:48
Diagnostic Test	Pass	08:03:15
Air Blank	0.000	08:03:43
0.100 Cal Check	0.097	08:04:02
Air Blank	0.000	08:04:32
Subject Test	0.250*	08:07:40
Air Blank	0.000	08:08:12
Air Blank	0.000	08:08:40
Subject Test	XXX**	08:09:06
Air Blank	0.000	08:09:40

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

03/11/2010

08:09:58

by


Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:10:25
Cal Check	0.096	08:10:44
Air Blank	0.000	08:11:14

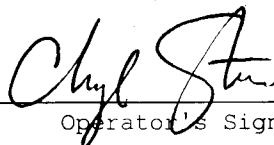

Operator's Signature

677
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
3-11-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 03/11/10 TIME 07:55
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

03/11/2010

07:55:24

154


Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:55:51
Cal Check	0.098	07:56:10
Air Blank	0.000	07:56:40



Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 1 SC

INTOXILYZER SERVICE RECORD

By
3.16.2010

Date: 2/26/10

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001165

Reason For Service: Breath hose was missing the part where a mouth piece attaches.

Date Out of Service: N/A

Calibration Check completed prior to removal from service

Date: N/A

SQAP completed prior to removal from service

Date: N/A

Corrective Action:

☐ Sent to CMI. Inc.

Replaced breath hose.