



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5-6-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/06/10 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 04/08/10 and 05/06/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

05/06/10


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6th, day of May, 2010.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

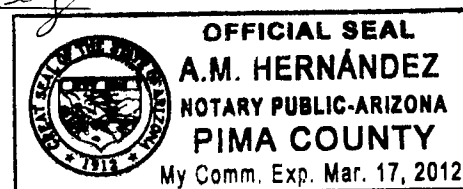

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
5.6.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 05/06/10 TIME 07:39
INTOXILYZER SERIAL # 80-001165 LOCATION TEAM 1

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

287
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

05/06/2010

07:41:20

BY

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:42:10
Diagnostic Test	Pass	07:42:37
Air Blank	0.000	07:43:05
0.100 Cal Check	0.099	07:43:23
Air Blank	0.000	07:43:54
Subject Test	0.000	07:44:25
Air Blank	ABT*	07:44:32
Air Blank	0.000	07:45:01

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

05/06/2010

07:45:09

TS4

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:45:40
Air Blank	IPS**	07:45:48

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

05/06/2010

07:45:59

134

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:46:49
Diagnostic Test	Pass	07:47:16
Air Blank	0.000	07:47:43
0.100 Cal Check	0.098	07:48:02
Air Blank	0.000	07:48:33
Subject Test	0.160*	07:51:41
Air Blank	0.000	07:52:14
Air Blank	0.000	07:52:41
Subject Test	XXX**	07:53:09
Air Blank	0.000	07:53:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

05/06/2010

07:53:58

TS4

Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:54:25
Cal Check	0.097	07:54:44
Air Blank	0.000	07:55:14



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

B4
5-6-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 05/06/10 TIME 07:39
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

05/06/2010

07:39:23

787
CS
154

Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:39:51
Cal Check	0.098	07:40:09
Air Blank	0.000	07:40:40



Operator's Signature