



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
7-6-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/01/10 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 06/03/10 and 07/01/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

07/01/10

ADPS QAS Permit # 33408


SUBSCRIBED AND SWORN to before me this 1, day of

July, 20 10.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

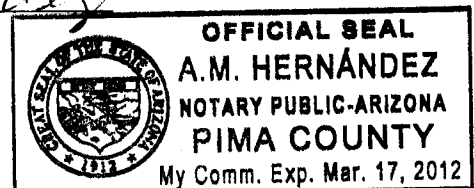

177
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
7-6-2010

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department

DATE 07/01/10 TIME 08:54

INTOXILYZER SERIAL # 80-001165 LOCATION TEAM 1

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

227
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
07/01/2010

08:56:08

TS4

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:56:59
Diagnostic Test	Pass	08:57:27
Air Blank	0.000	08:57:54
0.100 Cal Check	0.098	08:58:13
Air Blank	0.000	08:58:43
Subject Test	0.000	08:59:06
Air Blank	ABT*	08:59:13
Air Blank	0.000	08:59:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

387
CS

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
07/01/2010

08:59:45

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:00:18
Air Blank	IPS**	09:00:25

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
07/01/2010

09:00:33

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:01:22
Diagnostic Test	Pass	09:01:50
Air Blank	0.000	09:02:17
0.100 Cal Check	0.099	09:02:36
Air Blank	0.000	09:03:06
Subject Test	0.143*	09:06:15
Air Blank	0.000	09:06:47
Air Blank	0.000	09:07:14
Subject Test	XXX**	09:07:39
Air Blank	0.000	09:08:13

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
CS

154

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/01/2010

09:08:29

Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:08:57
Cal Check	0.097	09:09:15
Air Blank	0.000	09:09:46


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/01/10 TIME 08:54
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

797
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/01/2010

08:54:07

54

Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:54:34
Cal Check	0.098	08:54:53
Air Blank	0.000	08:55:22



Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 07/01/10
QAS: Cheryl STRAUB
Intoxilyzer Serial Number: 80-001165
New Lot Number: 719919
New Expiration Date: 02/02/2012

- 172 CS
154
7-6-2010
- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
 - (✓) 2. Remove the old gas tank and install the new one.
 - (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
 - (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
 - (✓) 5. Hit the "ESC" key twice.
 - (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
 - (✓) 7. Air blank completed.
 - (✓) 8. Calibration check completed. Test results 0. 098 AC.
 - (✓) 9. Air blank completed.
 - (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

277
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

07/01/2010

08:51:49

104


Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:52:16
Cal Check	0.098	08:52:35
Air Blank	0.000	08:53:05

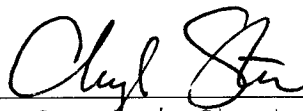

Operator's Signature