



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By
4.18.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 04/07/11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 03/10/2011 and 04/07/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

04/18/2011


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 18th, day of April, 20 11.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

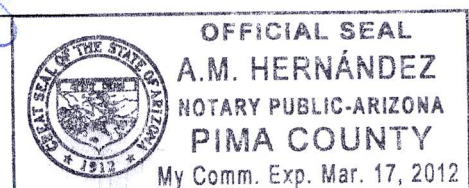

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

137
CSB4
4.13.2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 04/07/2011 TIME 07:39
INTOXILYZER SERIAL # 80-001165 LOCATION TEAM 1

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

277
CS

BH

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

04/07/2011

07:40:37

Standard Lot#: 719919

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:41:28
Diagnostic Test	Pass	07:41:56
Air Blank	0.000	07:42:23
0.100 Cal Check	0.099	07:42:42
Air Blank	0.000	07:43:12
Subject Test	0.000	07:43:33
Air Blank	ABT*	07:43:40
Air Blank	0.000	07:44:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
04/07/2011 07:44:13

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:44:46
Air Blank	IPS**	07:44:54

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001165
04/07/2011 07:45:09

Standard Lot#: 719919
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:45:59
Diagnostic Test	Pass	07:46:26
Air Blank	0.000	07:46:54
0.100 Cal Check	0.098	07:47:12
Air Blank	0.000	07:47:43
Subject Test	0.137*	07:50:51
Air Blank	0.000	07:51:23
Air Blank	0.000	07:51:51
Subject Test	XXX**	07:52:17
Air Blank	0.000	07:52:51

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
68

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

04/07/2011

07:53:09

154


Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:53:36
Cal Check	0.098	07:53:55
Air Blank	0.000	07:54:25


Operator's Signature

697
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

154
4.13.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 04/07/2011 TIME 07:39
Intoxilyzer Serial # 80-001165 Location TEAM 1

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

04/07/2011

07:39:15

104


Standard Lot#: 719919

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:39:42
Cal Check	0.098	07:40:01
Air Blank	0.000	07:40:31

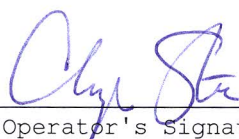

Operator's Signature