



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
7-6-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/30/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 06/02/11 and 06/30/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

07/6/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6, day of

July, 20 11.

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

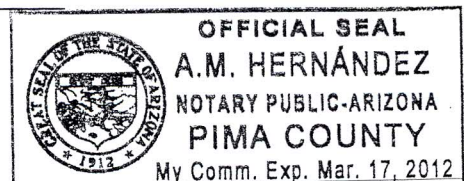


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

187  
CS  
B4  
7.6.2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 06/30/2011 TIME 07:16  
INTOXILYZER SERIAL # 80-001165 LOCATION TEAM I

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheryl Straub

287  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

06/30/2011

07:18:01

PH

Standard Lot#: 719919

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:18:49
Diagnostic Test	Pass	07:19:17
Air Blank	0.000	07:19:44
0.100 Cal Check	0.097	07:20:03
Air Blank	0.000	07:20:33
Subject Test	0.000	07:20:59
Air Blank	ABT*	07:21:05
Air Blank	0.000	07:21:33

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

387  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
06/30/2011

07:21:37

BY

---

Standard Lot#: 719919  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	07:22:08
Air Blank	IPS**	07:22:15

\*RFI Detect

\*\*Improper Sample

---

Not a Successfully Completed Test Sequence



487  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

06/30/2011

07:22:20

BH

Standard Lot#: 719919

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:23:10
Diagnostic Test	Pass	07:23:37
Air Blank	0.000	07:24:05
0.100 Cal Check	0.098	07:24:23
Air Blank	0.000	07:24:54
Subject Test	0.078*	07:28:02
Air Blank	0.000	07:28:34
Air Blank	0.000	07:29:02
Subject Test	XXX**	07:29:28
Air Blank	0.000	07:30:02

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

587  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

06/30/2011

07:30:25

BS4


-----  
Standard Lot#: 719919

Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:30:52
Cal Check	0.097	07:31:11
Air Blank	0.000	07:31:41

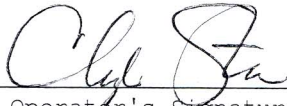  
\_\_\_\_\_  
Operator's Signature

687  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

06/30/2011

07:16:34

-----  
Standard Lot#: 719919

Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Bh

Test	g/210L	Time
Air Blank	0.000	07:17:01
Cal Check	0.097	07:17:20
Air Blank	0.000	07:17:50

  
Operator's Signature

787  
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

154  
7.6.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 06/30/11 TIME 07:16  
Intoxilyzer Serial # 80-001165 Location TEAM I

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)