



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
9.27-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.23-11 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 8.25-11 and 9.23-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

9.23-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of September, 2011.

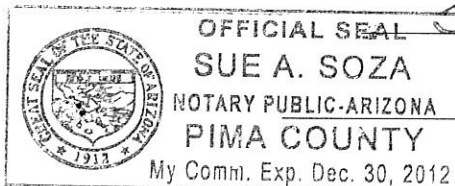

Sue A Soza

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1067

cu

By  
9.27-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 9.23-11 TIME 0816  
INTOXILYZER SERIAL # 80-001165 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001165

09/23/2011

08:11:56

Standard Lot#: 04911100A1

Last Changed By: C. STRAUB

#50772

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:12:39
Diagnostic Test	Pass	08:13:06
Air Blank	0.000	08:13:34
0.100 Cal Check	0.096	08:13:52
Air Blank	0.000	08:14:23
Subject Test	0.000	08:14:41
Air Blank	ABT*	08:14:47
Air Blank	0.000	08:15:16

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
09/23/2011

08:15:20

Standard Lot#: 04911100A1  
Last Changed By: C. STRAUB #50772  
QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:15:45
Air Blank	IPS**	08:15:54

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001165  
09/23/2011

08:16:01

467  
Cu

134

Standard Lot#: 04911100A1  
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:16:47
Diagnostic Test	Pass	08:17:15
Air Blank	0.000	08:17:42
0.100 Cal Check	0.096	08:18:01
Air Blank	0.000	08:18:31
Subject Test	0.408*	08:21:40
Air Blank	0.000	08:22:13
Air Blank	0.000	08:22:41
Subject Test	XXX**	08:23:01
Air Blank	0.000	08:23:35

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

09/23/2011

08:24:04

-----  
Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:24:31
Cal Check	0.094	08:24:50
Air Blank	0.000	08:25:21

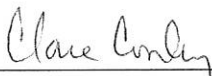
  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 9-23-11 TIME 0810  
Intoxilyzer Serial # 80-001165 Location ODS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

09/23/2011

08:10:29


-----  
Standard Lot#: 04911100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:10:57
Cal Check	0.095	08:11:15
Air Blank	0.000	08:11:46

  
\_\_\_\_\_  
Operator's Signature