



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TJ
11-19-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/8/12 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 10/11/12 and 11/8/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


11/14/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 12th, day of November, 2012.

Edma Escarcega

Notary Public

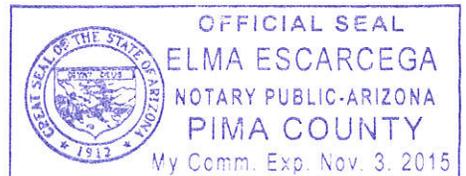

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY
11.19.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/8/12 TIME 1400

INTOXILYZER SERIAL # 80-001165 LOCATION Vehicle #945

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

BA

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/08/2012 14:02:25

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP,SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:03:22
Diagnostic Test	Pass	14:03:59
Air Blank	0.000	14:04:27
0.100 Cal Check	0.101	14:04:48
Air Blank	0.000	14:05:17
Subject Test	0.000	14:05:55
Air Blank	0.000	14:06:22
Wait	ABT*	
Air Blank	0.000	14:06:57

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7sc
BSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/08/2012 14:07:20

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:08:13
Diagnostic Test	Pass	14:08:50
Air Blank	0.000	14:09:18
0.100 Cal Check	0.102	14:09:38
Air Blank	RFI*	14:09:59
Air Blank	IPS**	14:10:07

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7_{sc}

BA

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/08/2012 14:10:16

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:11:05
Diagnostic Test	Pass	14:11:42
Air Blank	0.000	14:12:10
0.100 Cal Check	0.102	14:12:30
Air Blank	0.000	14:13:00
Subject Test	0.137*	14:16:07
Air Blank	0.000	14:16:39
Air Blank	0.000	14:17:07
Subject Test	XXX**	14:17:23
Air Blank	0.000	14:17:57

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7_{sc}

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/08/2012 14:18:16

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:18:44
0.100 Cal Check	0.101	14:19:05
Air Blank	0.000	14:19:35

Preventative Maintenance
Performed -- Timer Reset

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TSK
11-19-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/8/12 TIME 1400
Intoxilyzer Serial # 80-001165 Location Vehicle # 945

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
11/08/2012 14:00:49

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:01:18
0.100 Cal Check	0.102	14:01:38
Air Blank	0.000	14:02:08

Preventative Maintenance
Performed -- Timer Reset