



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1-28-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/3/13 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 12/6/12 and 1/3/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/7/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 7th, day of January, 2013.

Elma Escarcega

Notary Public

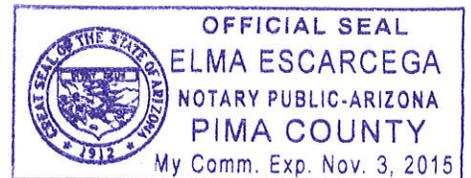

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc
B4
1.28.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/3/13 TIME 1048

INTOXILYZER SERIAL # 80-001165 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
01/03/2013 10:51:08

2 of 7 sc
Dh

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:00
Diagnostic Test	Pass	10:52:35
Air Blank	0.000	10:53:02
0.100 Cal Check	0.102	10:53:23
Air Blank	0.000	10:53:53
Subject Test	0.000	10:54:17
Air Blank	0.000	10:54:44
Wait	ABT*	
Air Blank	0.000	10:55:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
01/03/2013

10:59:02

3 of 7 sc
B4

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:59:53
Diagnostic Test	Pass	11:00:28
Air Blank	0.000	11:00:55
0.100 Cal Check	0.102	11:01:16
Air Blank	0.000	11:01:46
Subject Test	RFI*	11:02:20
Air Blank	IPS**	11:02:27

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
01/03/2013

11:02:35

4 of 7 sc
B4

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:26
Diagnostic Test	Pass	11:04:01
Air Blank	0.000	11:04:28
0.100 Cal Check	0.102	11:04:49
Air Blank	0.000	11:05:19
Subject Test	0.133*	11:08:26
Air Blank	0.000	11:08:58
Air Blank	0.000	11:09:26
Subject Test	XXX**	11:10:34
Air Blank	0.000	11:11:08

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

Core Version: 8105.48

01/03/2013

11:11:43

5 of 7 sc
B4

Standard Lot#: 13712100A2

Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:12:12
0.100 Cal Check	0.102	11:12:33
Air Blank	0.000	11:13:03

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
TS4
1.28.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/3/13 TIME 1048
Intoxilyzer Serial # 80-001165 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
01/03/2013

10:48:27

7 of 7 sc
124

Standard Lot#: 13712100A2
Last Changed By: W. HONOMICHL #42992

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:48:56
0.100 Cal Check	0.102	10:49:17
Air Blank	0.000	10:49:46

Preventative Maintenance
Performed -- Timer Reset