



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

BSy
9.3.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/15/14 Cal 10.100 AC Cal 2 N/A AC

Between the dates of 7/18/14 and 8/15/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 8/15/14.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


8/27/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27th, day of August, 2014.

Elma Escarcega

Notary Public

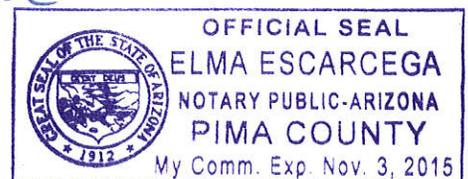

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 4
sc
TS4
28847
9.3.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/15/14 TIME 1419
Intoxilyzer Serial # 80-001165 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
08/15/2014 14:19:48

2 of 4 SL
BS

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:20:16
0.100 Cal Check	0.100	14:20:37
Air Blank	0.000	14:21:06

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
 CRIME LABORATORY
 INTOXILYZER SERVICE RECORD

3 of 4 sc
 TX
 9.3.2014

Date: 8/15/14

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001165

Reason For Service: "START TEST" button will not work to start subject test
Error message of "System Fault ID#: 011, Contact CMI" appeared.

Date Out of Service: 8/15/14

Calibration Check completed prior to removal from service

Date: 8/15/14

SQAP completed prior to removal from service

Date: N/A

Corrective Action:

Sent to CMI, Inc.



Service Evaluation Form

This form MUST be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

4 of 4 SL
TSy
288847
09 9.3.2014

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
Name Terry Gallegos Phone: (520) 837-7677
Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-001165 Instrument Model: I8000
appears after PM's reset (with arrow pointing to Serial Number)

4. Detailed Description of Problem:
After successful calibration check, message "System Fault ID#: 011, Contact CMI" appeared. Cannot start a subject test - START TEST button will not work, but does do initialization when instrument turned on. Also keyboard broken.

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator
Name (Please Print) Title
Terry Gallegos 8.22.2014
Signature Date

Ship to:
CMI, Inc.
Attn: Service Dept.
316 East Ninth Street
Owensboro, KY 42303

No, please send estimate before repairs are made.
Note: An estimate will be faxed before performing any repairs and may cause delays in service.
An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.