



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten initials and date:* 7/2/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.18.2015 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 5.26.2015 and 6.18.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

*Handwritten signature of Terry S. Gallegos*

Terry S. Gallegos


6.24.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24th, day of June, 2015.

*Handwritten signature of Elma Escarcega*

Notary Public

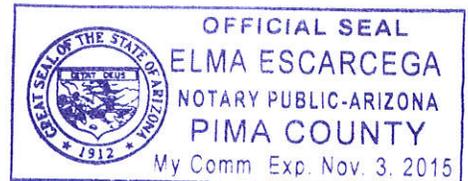

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8  
Tch  
Wb  
7/2/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 6-18-2015 TIME 1410

INTOXILYZER SERIAL # 80-001165 LOCATION DAI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015 14:14:09

2 of  
154  
JD

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:14:56
Diagnostic Test	Pass	14:15:31
Air Blank	0.000	14:15:59
0.100 Cal Check	0.097	14:16:20
Air Blank	0.000	14:16:49
Subject Test	0.000	14:17:04
Air Blank	0.000	14:17:31
Wait	ABT*	
Air Blank	0.000	14:18:13

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015 14:18:43

3 of 8  
754  
JJA

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:20:23
Diagnostic Test	Pass	14:20:58
Air Blank	0.000	14:21:25
0.100 Cal Check	0.097	14:21:46
Air Blank	RFI*	14:21:53
Air Blank	IPS**	14:22:08

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015

14:23:03

40f8  
TS4  
JL

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:23:48
Diagnostic Test	Pass	14:24:23
Air Blank	0.000	14:24:50
0.100 Cal Check	0.097	14:25:11
Air Blank	0.000	14:25:41
Subject Test	0.171*	14:28:48
Air Blank	0.000	14:29:20
Air Blank	0.000	14:29:48
Subject Test	0.086**	14:30:02
Air Blank	0.000	14:30:34
Wait		
Air Blank	0.000	14:35:30
Subject Test	0.000**	14:36:06
Air Blank	0.000	14:36:34

\*Deficient Sample

\*\*No .020 agreement

Subject test #2 was for Invalid sample test but I ran out of breath during the attempt. I waited 5min then did another test, followed by the Invalid Sample.  
TS4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015 14:36:59

*Sof8  
TSg  
JL*

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:37:49
Diagnostic Test	Pass	14:38:24
Air Blank	0.000	14:38:51
0.100 Cal Check	0.097	14:39:12
Air Blank	0.000	14:39:41
Subject Test	0.219	14:40:25
Air Blank	0.000	14:40:58
Wait		
Air Blank	0.000	14:45:52
Subject Test	0.000	14:48:32
Air Blank	0.000	14:48:59
Wait		
Air Blank	0.000	14:54:00
Subject Test	XXX*	14:54:46
Air Blank	0.000	14:55:20

\*Invalid Sample  
Begin new deprivation period

*On Subject test #1, the Invalid sample attempt, I ran out of breath during the attempt. I waited 5min, then followed with the Invalid Sample. TSg*  
**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015 14:55:34

60f8  
TSy  
Dd

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:56:03
0.100 Cal Check	0.097	14:56:23
Air Blank	0.000	14:56:53

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

70f8  
704  
JLS  
7/2/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6-18-2015 TIME 1410  
Intoxilyzer Serial # 80-001165 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

80f8  
D4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
06/18/2015 14:10:18

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:12:01
0.100 Cal Check	0.097	14:12:22
Air Blank	0.000	14:12:51

Preventative Maintenance  
Performed -- Timer Reset