



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

AS
 12-18-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12.3.2015 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 11.5.2015 and 12.3.2015 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Terry S. Gallegos
 Terry S. Gallegos

12-17-2015
 ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2015.

Edmo Escarcega
 Notary Public

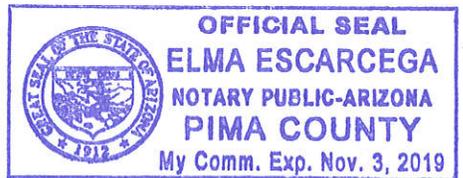

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 8
TS4
M
12-18-15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-3-2015 TIME 1254

INTOXILYZER SERIAL # 80-001165 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 800 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2018
TJ
AD

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015 13:05:06

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:05:58 |
| Diagnostic Test | Pass | 13:06:33 |
| Air Blank | 0.000 | 13:07:00 |
| 0.100 Cal Check | 0.100 | 13:07:21 |
| Air Blank | 0.000 | 13:07:50 |
| Subject Test | 0.000 | 13:08:24 |
| Air Blank | 0.000 | 13:08:51 |
| Wait | ABT* | |
| Air Blank | 0.000 | 13:09:32 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015 13:09:38

30f8
TS4
[Signature]

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:11:26 |
| Diagnostic Test | Pass | 13:12:01 |
| Air Blank | 0.000 | 13:12:28 |
| 0.100 Cal Check | RFI* | 13:12:42 |
| Air Blank | IPS** | 13:12:50 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015

13:13:32

4 of 8
TJh
[Signature]

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:14:23 |
| Diagnostic Test | Pass | 13:14:59 |
| Air Blank | 0.000 | 13:15:26 |
| 0.100 Cal Check | 0.099 | 13:15:47 |
| Air Blank | 0.000 | 13:16:16 |
| Subject Test | 0.125* | 13:19:24 |
| Air Blank | 0.000 | 13:19:56 |
| Air Blank | 0.000 | 13:20:23 |
| Subject Test | XXX** | 13:20:41 |
| Air Blank | 0.000 | 13:21:14 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015 13:21:37

5078
TS4
AD

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 13:21:59 |
| Air Blank | RFI* | 13:22:07 |

*RFI Detect

unintentional RFI TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015 13:22:20

60F8
TS4
Ad

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 13:22:49 |
| 0.100 Cal Check | 0.099 | 13:23:10 |
| Air Blank | 0.000 | 13:23:40 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

70f8
DB
Alb
12-18-15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12.3.2015 TIME 1254
Intoxilyzer Serial # 80-001165 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001165
Core Version: 8105.48
12/03/2015 12:54:36

8058
TJ
DD

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:55:06 |
| 0.100 Cal Check | 0.099 | 12:55:27 |
| Air Blank | 0.000 | 12:55:56 |
| Air Blank | 0.000 | 12:58:25 |

*Input Timed Out