



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten:* JLS  
5-26-16

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001165 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-19-2016 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 4-21-2016 and 5-19-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

*Terry S. Gallegos*

Terry S. Gallegos


5-24-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24th, day of May, 2016.

*Elma Escarcega*

Notary Public

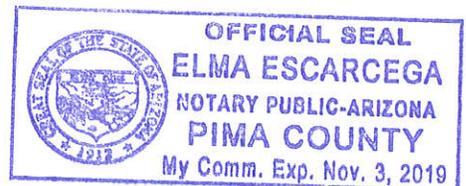

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1057  
TS4  
DLS  
5-26-16

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5-19-2016 TIME 1147

INTOXILYZER SERIAL # 80-001165 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
05/19/2016 11:50:04

2087  
Dh  
Ad

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:50:49
Diagnostic Test	Pass	11:51:24
Air Blank	0.000	11:51:51
0.100 Cal Check	0.098	11:52:12
Air Blank	0.000	11:52:42
Subject Test	0.000	11:52:56
Air Blank	0.000	11:53:23
Wait	ABT*	
Air Blank	0.000	11:53:58

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
05/19/2016

11:54:03

30F7  
TS4  
JAD

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:55:29
Diagnostic Test	Pass	11:56:04
Air Blank	0.000	11:56:32
0.100 Cal Check	RFI*	11:56:47
Air Blank	IPS**	11:56:55

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
05/19/2016

11:57:34

4 of 7  
DL  
JL

---

Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:58:24
Diagnostic Test	Pass	11:58:59
Air Blank	0.000	11:59:26
0.100 Cal Check	0.098	11:59:47
Air Blank	0.000	12:00:17
Subject Test	0.243*	12:03:24
Air Blank	0.000	12:03:56
Air Blank	0.000	12:04:24
Subject Test	XXX**	12:05:22
Air Blank	0.000	12:05:55

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001165

Core Version: 8105.48

05/19/2016

12:07:43

507  
TS4  
[Signature]

-----  
Standard Lot#: 15615100A1

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD  
-----

Operator: T. GALLEGOS #28847

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:08:13
0.100 Cal Check	0.098	12:08:34
Air Blank	0.000	12:09:03

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
TS4  
DDE  
5-26-16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 5.19.2016 TIME 1147  
Intoxilyzer Serial # 80-001165 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001165  
Core Version: 8105.48  
05/19/2016 11:47:27

7057  
TS4  
JLD

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:47:56
0.100 Cal Check	0.098	11:48:17
Air Blank	0.000	11:48:46

Preventative Maintenance  
Performed -- Timer Reset