



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

jk

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-28-2006 Cal 1 0.097 AC Cal 2 0.094 AC

Between the dates of 12-1-2006 and 12-28-2006 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos
Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

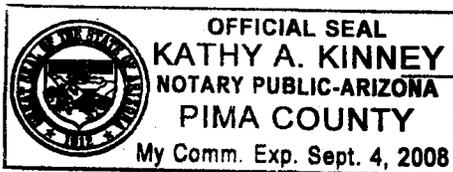

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1078
TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

ju

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28047 AGENCY Tucson Police Department

DATE 12-28-2006 TIME 1742

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001199
12/28/2006

ju

17:43:46

2078
754

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:44:37
Diagnostic Test	Pass	17:45:05
Air Blank	0.000	17:45:32
0.100 Cal Check	0.096	17:45:51
Air Blank	0.000	17:46:20
Subject Test	0.000	17:46:37
Air Blank	0.000	17:47:04
Five Minute Wait	ABT*	
Air Blank	0.000	17:47:53

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001199
12/28/2006

je

17:48:20

3 of 8
D4

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:49:09
Diagnostic Test	Pass	17:49:36
Air Blank	IPS*	17:49:53
Air Blank	0.000	17:50:21

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001199
12/28/2006 17:50:46

4058
TS4

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:51:37
Diagnostic Test	Pass	17:52:04
Air Blank	0.000	17:52:31
0.100 Cal Check	RFI*	17:52:46
Air Blank	0.000	17:53:16

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001199
12/28/2006

je

17:53:46

50FB
TB4

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:54:47
Diagnostic Test	Pass	17:55:15
Air Blank	0.000	17:55:42
0.100 Cal Check	0.097	17:56:01
Air Blank	0.000	17:56:30
Subject Test	0.059*	17:59:39
Air Blank	0.000	18:00:11
Air Blank	0.000	18:00:38
Subject Test	XXX**	18:02:59
Air Blank	0.000	18:03:34

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001199
12/28/2006 18:04:06

je

6 of 8
TD4

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	18:04:34
Cal Check	0.094	18:04:52
Air Blank	0.000	18:05:22


Operator's signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

70f8
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

jk

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 2004 AGENCY Tucson Police Department
DATE 12-28-2006 TIME 1742
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001199
12/28/2006 17:42:19

gk

8058
754

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:42:46
Cal Check	0.097	17:43:05
Air Blank	0.000	17:43:34


Operator's Signature