



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

7m


I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-26-2007 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 1-26-2007 and 2-26-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

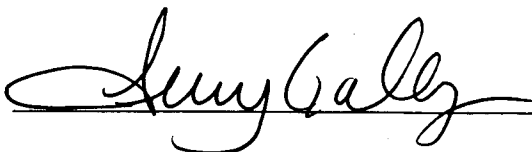

Terry S. Gallegos

March 12, 2007


ADPS QAS Permit # 33169



SUBSCRIBED AND SWORN to before me this 12, day of

OFFICIAL SEAL
KATHY A. KINNEY
NOTARY PUBLIC-ARIZONA
PIMA COUNTY
My Comm. Exp. Sept. 4, 2008

March, 2007.

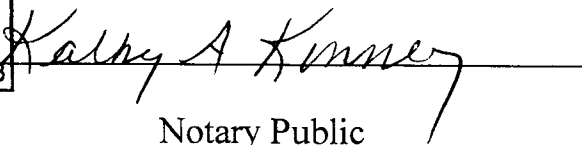
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1078
PSY ms

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 2-26-2007 TIME 1038

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/26/2007

10:47:29

2 of 8
Dh m

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:48:21
Diagnostic Test	Pass	10:48:49
Air Blank	0.000	10:49:16
0.100 Cal Check	0.099	10:49:35
Air Blank	0.000	10:50:04
Subject Test	0.000	10:50:34
Air Blank	ABT*	10:50:50
Air Blank	0.000	10:51:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/26/2007

10:51:47

30f 8
Bh mm

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:41
Diagnostic Test	Pass	10:53:09
Air Blank	RFI*	10:53:18
Air Blank	IPS**	10:53:27

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

02/26/2007

10:53:34

40f8
B4 m

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS

#28847

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:54:29
Diagnostic Test	Fail*	10:54:33
Air Blank	0.000	10:55:01

* Analytical Stability Test

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
02/26/2007

10:55:59

508
By MM

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:56:49
Diagnostic Test	Pass	10:57:16
Air Blank	0.000	10:57:43
0.100 Cal Check	0.098	10:58:02
Air Blank	0.000	10:58:31
Subject Test	0.123*	11:01:40
Air Blank	0.000	11:02:11
Air Blank	0.000	11:02:39
Subject Test	XXX**	11:03:09
Air Blank	0.000	11:03:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

02/26/2007

11:04:15

60f 8
Tsh ms

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:04:43
Cal Check	0.098	11:05:02
Air Blank	0.000	11:05:31


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8
Dy me

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28049 AGENCY Tucson Police Department
DATE 2-26-2007 TIME 1038
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

02/26/2007

10:37:44

Bof 8
Bhms

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	10:38:11
Cal Check	0.100	10:38:30
Air Blank	0.000	10:38:59


Operator's Signature