



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

*M*  
8.22-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.19.2007 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 6.21.2007 and 7.19.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.


*Terry S. Gallegos*

Terry S. Gallegos

8.15.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15, day of August, 2007.

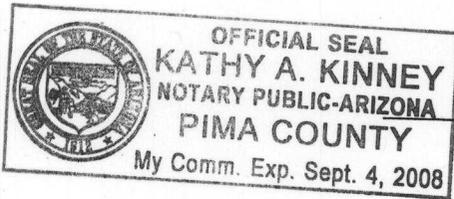

*Kathy A Kinney*  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

CU 8-22-07  
1 of 7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7-19-2007 TIME 1827

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
07/19/2007 18:28:56

du  
2 of 7  
BY

---

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	18:29:47
Diagnostic Test	Pass	18:30:14
Air Blank	0.000	18:30:42
0.100 Cal Check	0.098	18:31:00
Air Blank	0.000	18:31:30
Subject Test	0.000	18:31:59
Air Blank	0.000	18:32:27
Five Minute Wait	ABT*	
Air Blank	0.000	18:33:03

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
07/19/2007 18:33:39

cu  
3077  
T34

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:34:30
Diagnostic Test	Pass	18:34:58
Air Blank	0.000	18:35:25
0.100 Cal Check	0.096	18:35:44
Air Blank	RFI*	18:35:53
Air Blank	IPS**	18:36:05

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
07/19/2007 18:36:23

cu  
4077  
154

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:38:32
Diagnostic Test	Pass	18:39:00
Air Blank	0.000	18:39:27
0.100 Cal Check	0.097	18:39:46
Air Blank	0.000	18:40:15
Subject Test	0.133*	18:43:24
Air Blank	0.000	18:43:55
Air Blank	0.000	18:44:22
Subject Test	XXX**	18:44:45
Air Blank	0.000	18:45:18

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
07/19/2007 18:45:54

cu  
5087  
TS4

-----  
Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:46:22
Cal Check	0.096	18:46:40
Air Blank	0.000	18:47:10

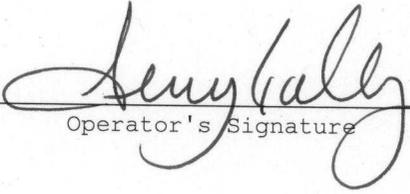
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)


ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

du 8-22-07  
6077  
TS4

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.19.2007 TIME 1827  
Intoxilyzer Serial # 80-001199 Location JUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

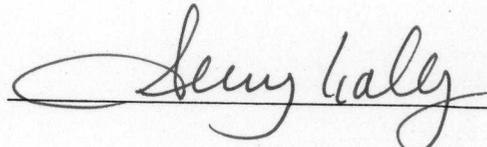SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
07/19/2007 18:27:35

OK  
7 of 7  
TS4

-----  
Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:28:02
Cal Check	0.096	18:28:21
Air Blank	0.000	18:28:50

  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

DM  
8.22-07

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 8.7.2007  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-001199  
New Lot Number: 627101G  
New Expiration Date: 10.03.2008

1 of 2  
TBH

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:  
 Enter Target Value: 0.100  
 Enter Control Lot #: From the tank label  
 Enter Exp Date: From the tank label (mm/dd/yyyy)  
 Edit info Y/N? Y or N  
 Update Std Info? Y ENTER  
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 097 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
08/07/2007

15:30:05

OK  
2 of 2  
B4

-----  
Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	15:30:33
Cal Check	0.097	15:30:52
Air Blank	0.000	15:31:21

  
-----  
Operator's Signature