



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/27/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 06/25/09 and 07/27/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Cheryl K. Straub

Cheryl K. Straub


07/28/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 28, day of July, 2009.

A.M. Hernández

Notary Public

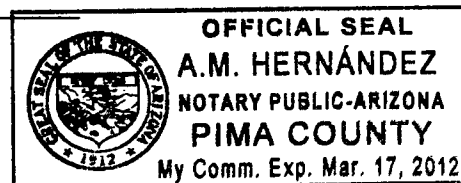

\*Includes documents specified in A.A.C. R13-10-104.A.5

197  
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDUREBY  
7-28-09QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/27/09 TIME 13:52  
INTOXILYZER SERIAL # 80-001199 LOCATION CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- 
- ( ✓ ) 2. Date check.

## OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0.
- 000
- AC.
- 
- ( ✓ ) 2. Error recognition logic system functioning.
- 
- Not a Successfully Completed Test Sequence printed
- 
- ( ✓ ) 3. Proper sample recognition system.
- 
- Not a Successfully Completed Test Sequence printed
- 
- Deficient sample printed.
- 
- ( ✓ ) 4. Standard Calibration Check standard 0.
- 100
- AC
- 
- Result: 0.
- 099
- AC

Instrument is operating properly and accurately. YES ☒ NO ☐

## COMMENTS

---

---

---

---

SIGNATURE

Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

07/27/2009

14:01:28

104

Standard Lot#: 801002I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:02:32
Diagnostic Test	Pass	14:03:00
Air Blank	0.000	14:03:27
0.100 Cal Check	0.101	14:03:46
Air Blank	0.000	14:04:15
Subject Test	0.000	14:04:40
Air Blank	ABT*	14:04:47
Air Blank	0.000	14:05:15

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

387 CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
07/27/2009

14:05:41

154

---

Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	14:06:23
Air Blank	IPS**	14:06:33

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

497  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

07/27/2009

14:06:56

154

Standard Lot#: 801002I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:08:05
Diagnostic Test	Pass	14:08:32
Air Blank	0.000	14:09:00
0.100 Cal Check	0.100	14:09:18
Air Blank	0.000	14:09:48
Subject Test	0.047*	14:12:56
Air Blank	0.000	14:13:27
Air Blank	0.000	14:13:55
Subject Test	XXX**	14:14:32
Air Blank	0.000	14:15:05

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

577  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

07/27/2009

14:15:45

154

-----  
Standard Lot#: 801002I

Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:16:12
Cal Check	0.099	14:16:31
Air Blank	0.000	14:17:00



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000Tdy  
7-28-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/27/09 TIME 13:52  
Intoxilyzer Serial # 80-001199 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

07/27/2009

13:52:30

104


-----  
Standard Lot#: 801002I

Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:52:58
Cal Check	0.100	13:53:16
Air Blank	0.000	13:53:46

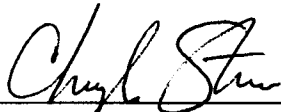

Operator's Signature