



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

134
8-31-2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/27/09 Cal 1 0.101 AC Cal 2 0.098 AC

Between the dates of 07/27/09 and 08/27/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

08/31/09


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 20 09.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

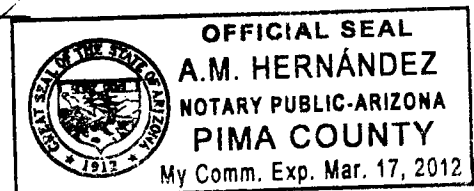

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
8.31.09

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 08/27/09 TIME 11:01
INTOXILYZER SERIAL # 80-001199 LOCATION CRIME LAB

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

(☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

287 cz

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

08/27/2009

11:12:18

134

Standard Lot#: 801002I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:13:11
Diagnostic Test	Pass	11:13:39
Air Blank	0.000	11:14:06
0.100 Cal Check	0.100	11:14:25
Air Blank	0.000	11:14:54
Subject Test	0.000	11:15:14
Air Blank	ABT*	11:15:21
Air Blank	0.000	11:15:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
cs

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

08/27/2009

11:15:55

BM

Standard Lot#: 801002I

Last Changed By: S. CIONE

#50616

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:16:37
Air Blank	IPS**	11:16:44

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

08/27/2009

11:16:54

Standard Lot#: 801002I

Last Changed By: S. CIONE

#50616

BY

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TPD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:17:59
Diagnostic Test	Pass	11:18:27
Air Blank	0.000	11:18:54
0.100 Cal Check	0.100	11:19:13
Air Blank	0.000	11:19:42
Subject Test	0.147*	11:22:50
Air Blank	0.000	11:23:23
Air Blank	0.000	11:23:50
Subject Test	XXX**	11:24:20
Air Blank	0.000	11:24:53

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

08/27/2009

11:25:19

Standard Lot#: 801002I


Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772

TUCSON PD

BM

Test	g/210L	Time
Air Blank	0.000	11:25:47
Cal Check	0.098	11:26:06
Air Blank	0.000	11:26:35

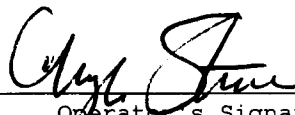

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

B4
8.21.09

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 08/27/09 TIME 11:01
Intoxilyzer Serial # 80-001199 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
c8

124


1

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
08/27/2009 11:01:37

Standard Lot#: 8010021
Last Changed By: S. CIONE #50616

OAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:02:05
Cal Check	0.101	11:02:24
Air Blank	0.000	11:02:53

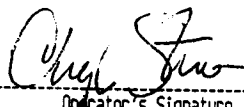
Operator's Signature