



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154
3-10-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-4-10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2-4-10 and 3-4-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


3-4-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 4th day of March, 2010.

Susan A. Hitchye

Notary Public

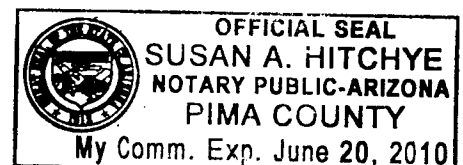

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
cu
B4
3.10.2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 3-4-10 TIME 0812

INTOXILYZER SERIAL # 80-001199 LOCATION ODW DUI HO

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None,

SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
03/04/2010

08:22:51

2 of 7
cu
B4

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:23:43
Diagnostic Test	Pass	08:24:10
Air Blank	0.000	08:24:37
0.100 Cal Check	0.099	08:24:56
Air Blank	0.000	08:25:25
Subject Test	0.000	08:25:50
Air Blank	ABT*	08:25:57
Air Blank	0.000	08:26:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
03/04/2010

08:26:30

367
en
/84

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:27:08
Air Blank	IPS**	08:27:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

03/04/2010

08:28:08

Standard Lot#: 642845

Last Changed By: J. OLIVARES

#36552

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:28:58
Diagnostic Test	Pass	08:29:25
Air Blank	0.000	08:29:52
0.100 Cal Check	0.099	08:30:11
Air Blank	0.000	08:30:41
Subject Test	0.165*	08:33:49
Air Blank	0.000	08:34:21
Air Blank	0.000	08:34:48
Subject Test	XXX**	08:35:35
Air Blank	0.000	08:36:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

03/04/2010

08:36:26

Standard Lot#: 642845

Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:36:54
Cal Check	0.099	08:37:12
Air Blank	0.000	08:37:42

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

657
cu

by
3.10.2010

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 3-4-10 TIME 0812
Intoxilyzer Serial # 80-001199 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

03/04/2010

08:21:12

Standard Lot#: 642845

Last Changed By: J. OLIVARES #36552

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:21:39
Cal Check	0.100	08:21:58
Air Blank	0.000	08:22:27

Clare Conley

Operator's Signature