



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/24/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 6/4/10 and 6/24/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

6/24/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24, day of June, 20 10.

Ant Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

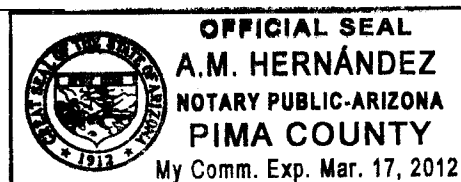

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
6-25-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 6/24/10 TIME 949INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Sarah Cione

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/24/2010

09:51:32

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:52:17
Diagnostic Test	Pass	09:52:44
Air Blank	0.000	09:53:12
0.100 Cal Check	0.100	09:53:31
Air Blank	0.000	09:54:00
Subject Test	0.000	09:54:35
Air Blank	0.000	09:55:02
Five Minute Wait	ABT*	
Air Blank	0.000	09:55:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/24/2010

09:56:59

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:58:10
Diagnostic Test	Pass	09:58:37
Air Blank	RFI*	09:59:01
Air Blank	IPS**	09:59:11

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/24/2010

09:59:27

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP,
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:00:28
Diagnostic Test	Pass	10:00:56
Air Blank	0.000	10:01:23
0.100 Cal Check	0.100	10:01:42
Air Blank	0.000	10:02:11
Subject Test	0.120*	10:05:20
Air Blank	0.000	10:05:51
Air Blank	0.000	10:06:19
Subject Test	XXX**	10:06:46
Air Blank	0.000	10:07:19

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

06/24/2010

10:07:38

Standard Lot#: 642845

Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:08:05
Cal Check	0.099	10:08:24
Air Blank	0.000	10:08:53

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
6-25-2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/24/10 TIME 949
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

06/24/2010

09:49:53

By

Standard Lot#: 642845

Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:20
Cal Check	0.100	09:50:39
Air Blank	0.000	09:51:08

Sarah Cione

Operator's Signature