



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
9-21-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/16/10 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8/19/10 and 9/16/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


9/20/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th day of September, 2010.

A.M. Hernandez

Notary Public

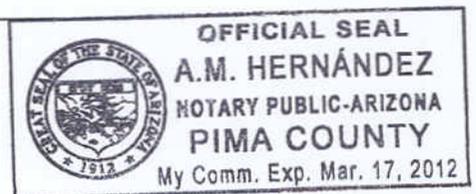

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

Tby
9.21.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/16/10 TIME 1123

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
09/16/2010

11:24:48

2 of 7sc

PH

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:25:45
Diagnostic Test	Pass	11:26:13
Air Blank	0.000	11:26:40
0.100 Cal Check	0.099	11:26:59
Air Blank	0.000	11:27:28
Subject Test	0.000	11:27:57
Air Blank	0.000	11:28:25
Five Minute Wait	ABT*	
Air Blank	0.000	11:29:20

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
09/16/2010 11:30:17

3 of 7 SC
154

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:31:03
Diagnostic Test	Pass	11:31:30
Air Blank	0.000	11:31:57
0.100 Cal Check	RFI*	11:32:11
Air Blank	IPS**	11:32:24

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
09/16/2010

11:32:42

4 of 7sc
AB4

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:33:26
Diagnostic Test	Pass	11:33:54
Air Blank	0.000	11:34:21
0.100 Cal Check	0.099	11:34:40
Air Blank	0.000	11:35:09
Subject Test	0.065*	11:38:17
Air Blank	0.000	11:38:49
Air Blank	0.000	11:39:16
Subject Test	XXX**	11:39:47
Air Blank	0.000	11:40:19

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
09/16/2010 11:40:37

5 of 7 sc

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

TSK

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:41:04
Cal Check	0.099	11:41:23
Air Blank	0.000	11:41:52

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 se

TDH
9.21.2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/16/10 TIME 1123
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
09/16/2010 11:23:01

7 of 75c

Standard Lot#: 642845
Last Changed By: J. OLIVARES #36552

184

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:23:28
Cal Check	0.099	11:23:47
Air Blank	0.000	11:24:17

Sarah Cione
Operator's Signature