



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY  
10-18-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/10 Cal 1 0.099 AC Cal 2 0.100 AC

Between the dates of 9/16/10 and 10/14/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


10/14/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14 day of October, 2010.

A.M. Hernandez

Notary Public

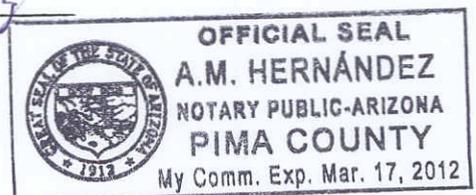

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY  
10-18-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/10 TIME 939

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 8 sc

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
10/14/2010 09:41:00

124

Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:44
Diagnostic Test	Pass	09:42:12
Air Blank	0.000	09:42:39
0.100 Cal Check	0.099	09:42:58
Air Blank	0.000	09:43:27
Subject Test	0.000	09:43:47
Air Blank	0.000	09:44:14
Five Minute Wait	ABT*	
Air Blank	0.000	09:44:50

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
10/14/2010

09:45:13

3 of 8 sc

BH

Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:45:54
Diagnostic Test	Pass	09:46:21
Air Blank	0.000	09:46:49
0.100 Cal Check	0.100	09:47:07
Air Blank	RFI*	09:47:34
Air Blank	0.000	09:48:02

\*RFI Detect

**Not a Successfully Completed Test Sequence**

4 of 8 SC

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
10/14/2010 09:48:55

BY

Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	09:49:43
Air Blank	0.000	09:50:12

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001199  
10/14/2010

09:50:24

5 of 8 sc

AS4

Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:51:06
Diagnostic Test	Pass	09:51:33
Air Blank	0.000	09:52:00
0.100 Cal Check	0.100	09:52:19
Air Blank	0.000	09:52:49
Subject Test	0.064*	09:55:57
Air Blank	0.000	09:56:28
Air Blank	0.000	09:56:56
Subject Test	XXX**	09:57:13
Air Blank	0.000	09:57:45

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
10/14/2010 09:58:05

6 of 8 sc

B4

-----  
Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:58:32
Cal Check	0.100	09:58:51
Air Blank	0.000	09:59:20

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

By  
10.18.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/14/10 TIME 939  
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
10/14/2010 09:39:16

8 of 8 sc

109

-----  
Standard Lot#: 642845  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:39:43
Cal Check	0.099	09:40:02
Air Blank	0.000	09:40:31

*Sarah Cione*

Operator's Signature