



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSH
6.27.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/23/11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 5/26/11 and 6/23/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

6/27/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of June, 20 11.

Am Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

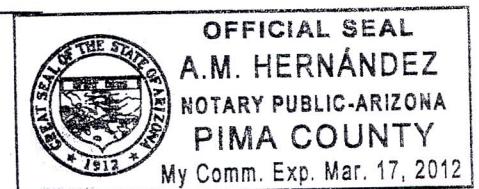

Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

124
6-27-2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 6/23/11 TIME 953INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001199
06/23/2011

09:55:08

Standard Lot#: 719919
Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:55:58
Diagnostic Test	Pass	09:56:25
Air Blank	0.000	09:56:53
0.100 Cal Check	0.100	09:57:11
Air Blank	0.000	09:57:41
Subject Test	0.000	09:58:36
Air Blank	0.000	09:59:03
Five Minute Wait	ABT*	
Air Blank	0.000	09:59:50

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

06/23/2011

10:02:59

Standard Lot#: 719919

Last Changed By: C. MORAND

#24219

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:03:23
Air Blank	IPS**	10:03:31

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

06/23/2011

10:06:25

Standard Lot#: 719919

Last Changed By: C. MORAND

#24219

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:07:22
Diagnostic Test	Pass	10:07:49
Air Blank	0.000	10:08:17
0.100 Cal Check	0.100	10:08:35
Air Blank	0.000	10:09:05
Subject Test	0.051*	10:12:13
Air Blank	0.000	10:12:44
Air Blank	0.000	10:13:11
Subject Test	XXX**	10:13:28
Air Blank	0.000	10:14:00

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

06/23/2011

10:14:31

Standard Lot#: 719919

Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:14:59
Cal. Check	0.099	10:15:17
Air Blank	0.000	10:15:47

Sarah Cione

Operator's Signature

5 of 9_{sc}

154

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/23/11 TIME 953
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
06/23/2011

09:53:06

Standard Lot#: 719919
Last Changed By: C. MORAND #24219

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:53:33
Cal Check	0.099	09:53:52
Air Blank	0.000	09:54:21



Operator's Signature

7 of 9_{sc}

134



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9_{sc}

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TDH
6-27-2011

Date: 6/23/11
QAS: Sarah Cione
Intoxilyzer Serial Number: 80-001199
New Lot Number: 20910100A1
New Expiration Date: 8/1/12

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START".
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing. "
Select "D" for dry control test. After selection is made press
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
06/23/2011 10:29:29

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:29:57
Cal Check	0.099	10:30:15
Air Blank	0.000	10:30:45

Sarah Cione
Operator's Signature