



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TJG
4.11.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/29/12 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 3/1/12 and 3/29/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 20 12.

Adina Escarcega

Notary Public

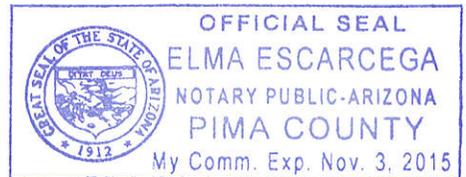

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
4.11.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 3/29/12 TIME 1349
INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/29/2012 13:51:43

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:52:28
Diagnostic Test	Pass	13:53:03
Air Blank	0.000	13:53:31
0.100 Cal Check	0.099	13:53:52
Air Blank	0.000	13:54:21
Subject Test	0.000	13:54:35
Air Blank	0.000	13:55:02
Wait	ABT*	
Air Blank	0.000	13:56:30

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7 sc

TS4

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-001199
 Core Version: 8105.48
 03/29/2012 13:56:42

Standard Lot#: 20910100A1
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:57:11
Air Blank	IPS**	13:57:18

*RFI Detect
 **Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/29/2012 14:03:01

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:03:58
Diagnostic Test	Pass	14:04:33
Air Blank	0.000	14:05:01
0.100 Cal Check	0.099	14:05:21
Air Blank	0.000	14:05:51
Subject Test	0.054*	14:08:58
Air Blank	0.000	14:09:30
Air Blank	0.000	14:09:57
Subject Test	XXX**	14:10:10
Air Blank	0.000	14:10:43

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/29/2012 14:13:00

5 of 7 sc

TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:13:29
0.100 Cal Check	0.098	14:13:50
Air Blank	0.000	14:14:20

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
4.11.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/29/12 TIME 1349
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/29/2012 13:49:16

TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:49:44
0.100 Cal Check	0.100	13:50:05
Air Blank	0.000	13:50:34

Preventative Maintenance
Performed -- Timer Reset