



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TH  
 4.16.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/29/13 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 2/28/13 and 3/29/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/8/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 8th, day of April, 2013.

Edma Escarcega

Notary Public

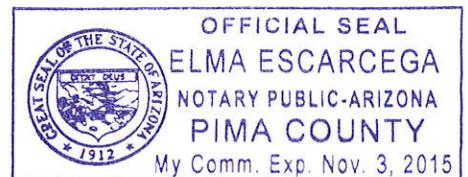

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7  
sc  
TB4  
4.16.2013

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/29/13 TIME 1527

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2077sc  
124

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/29/2013 15:28:54

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:29:36
Diagnostic Test	Pass	15:30:11
Air Blank	0.000	15:30:38
0.100 Cal Check	0.097	15:30:59
Air Blank	0.000	15:31:28
Subject Test	0.000	15:31:43
Air Blank	0.000	15:32:10
Wait	ABT*	
Air Blank	0.000	15:32:44

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7  
134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/29/2013 15:32:49

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON POLICE DEPT

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	15:34:05
Air Blank	IPS**	15:34:14

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 sc  
184

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/29/2013 15:34:24

---

Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	15:35:12
Diagnostic Test	Pass	15:35:47
Air Blank	0.000	15:36:15
0.100 Cal Check	0.099	15:36:36
Air Blank	0.000	15:37:05
Subject Test	0.161*	15:40:12
Air Blank	0.000	15:40:44
Air Blank	0.000	15:41:11
Subject Test	XXX**	15:41:27
Air Blank	0.000	15:42:01

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7  
TB

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/29/2013 15:42:59

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:43:28
0.100 Cal Check	0.098	15:43:48
Air Blank	0.000	15:44:18

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7  
TS4  
4.16.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/29/13 TIME 1527  
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 SC  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/29/2013 15:27:26

-----  
Standard Lot#: 23511100A2  
Last Changed By: J. HAWKINS #42401  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:27:54
0.100 Cal Check	0.099	15:28:15
Air Blank	0.000	15:28:44

Preventative Maintenance  
Performed -- Timer Reset