



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
6-5-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/23/13 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 4/25/13 and 5/23/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


5/30/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30th, day of May, 2013.

Elma Escarcega

Notary Public

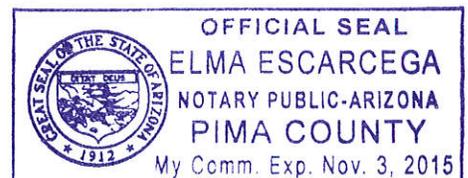

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
TS4
6-5-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 5/23/13 TIME 925
INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
05/23/2013

09:27:33

2 of 7 sc
Sh

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:44
Diagnostic Test	Pass	09:29:19
Air Blank	0.000	09:29:47
0.100 Cal Check	0.100	09:30:08
Air Blank	0.000	09:30:37
Subject Test	0.000	09:31:04
Air Blank	0.000	09:31:32
Wait	ABT*	
Air Blank	0.000	09:32:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
05/23/2013 09:32:13

3 of 7 sec
154

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:33:18
Diagnostic Test	Pass	09:33:53
Air Blank	0.000	09:34:21
0.100 Cal Check	0.099	09:34:41
Air Blank	0.000	09:35:11
Subject Test	RFI*	09:35:38
Air Blank	IPS**	09:35:46

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
05/23/2013

09:42:39

4 of 7
sc
10h

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:43:33
Diagnostic Test	Pass	09:44:08
Air Blank	0.000	09:44:36
0.100 Cal Check	0.098	09:44:56
Air Blank	0.000	09:45:26
Subject Test	0.176*	09:48:33
Air Blank	0.000	09:49:05
Air Blank	0.000	09:49:32
Subject Test	XXX**	09:50:00
Air Blank	0.000	09:50:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
05/23/2013 09:50:43

5 of 7
sc
TS4

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:51:11
0.100 Cal Check	0.098	09:51:32
Air Blank	0.000	09:52:02

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

607 SC
TS4
6-5-2013

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/23/13 TIME 925
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
05/23/2013 09:25:20

7 of 7 SC
TS4

Standard Lot#: 23511100A2
Last Changed By: J. HAWKINS #42401

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:25:48
0.100 Cal Check	0.099	09:26:09
Air Blank	0.000	09:26:38

Preventative Maintenance
Performed -- Timer Reset