



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSY
12.3.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/7/13 Cal 1 0.096 AC Cal 2 0.097 AC

Between the dates of 10/10/13 and 11/7/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


11/20/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th day of November, 2013.

Elma Escarcega

Notary Public

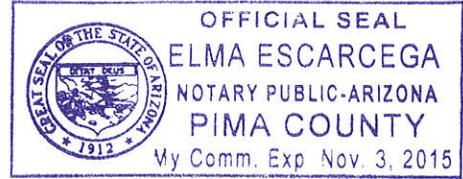

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
By
12-3-2013

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/7/13 TIME 1006

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/07/2013

10:09:48

2 of 7 sc
B4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:10:47
Diagnostic Test	Pass	10:11:22
Air Blank	0.000	10:11:49
0.100 Cal Check	0.096	10:12:10
Air Blank	0.000	10:12:40
Subject Test	0.000	10:13:17
Air Blank	0.000	10:13:44
Wait	ABT*	
Air Blank	0.000	10:15:30

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/07/2013

10:15:51

3 of 7 sc
B4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:16:45
Diagnostic Test	Pass	10:17:20
Air Blank	0.000	10:17:47
0.100 Cal Check	0.097	10:18:08
Air Blank	0.000	10:18:37
Subject Test	RFI*	10:18:56
Air Blank	IPS**	10:19:03

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/07/2013

10:23:02

4 of 7
SC
B4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:23:56
Diagnostic Test	Pass	10:24:31
Air Blank	0.000	10:24:59
0.100 Cal Check	0.098	10:25:20
Air Blank	0.000	10:25:49
Subject Test	0.057*	10:28:56
Air Blank	0.000	10:29:28
Air Blank	0.000	10:29:55
Subject Test	XXX**	10:31:14
Air Blank	0.000	10:31:47

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/07/2013 10:32:11

5 of 7 sc
DL

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:32:40
0.100 Cal Check	0.097	10:33:01
Air Blank	0.000	10:33:30

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TS4
12-3-2013

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/7/13 TIME 1006
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/07/2013 10:06:47

7 of 7 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:07:15
0.100 Cal Check	0.096	10:07:36
Air Blank	0.000	10:08:05

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2
sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

PK
12-3-2013

Date: 11/2/13
QAS: SULLIVAN QAS # 3327Z
Intoxilyzer Serial Number: 80-001199
New Lot Number: 30412100A1
New Expiration Date: 12/1/14

- () 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 2. Remove the old gas tank and install the new one.
- () 3. Hit the "ESC" key twice and swipe the QAS card.
- () 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- () 5. Hit the "ESC" key twice.
- () 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- () 7. Air blank completed.
- () 8. Calibration check completed. Test results 0. 095 AC.
- () 9. Air blank completed.
- () 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE 

2 of 2sc
TS4

SULLIVAN
SULLIVAN # 33272

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
11/02/2013 20:43:52

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. SULLIVAN #34971
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	20:44:26
0.100 Cal Check	0.095	20:44:47
Air Blank	0.000	20:45:16

Preventative Maintenance NOT
Performed -- Timer NOT Reset