



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*TSy*  
 2-12-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/30/14 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 1/2/14 and 1/30/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/11/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11<sup>th</sup>, day of February, 2014.

Elma Escarcega

Notary Public

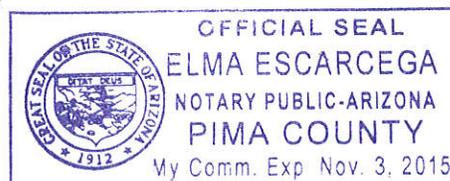

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc  
Bh  
2-12-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/30/14 TIME 839

INTOXILYZER SERIAL # 80-001199 LOCATION DUT HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
01/30/2014

08:41:29

2 of 7 sc  
TSK

---

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:42:16
Diagnostic Test	Pass	08:42:51
Air Blank	0.000	08:43:18
0.100 Cal Check	0.096	08:43:39
Air Blank	0.000	08:44:08
Subject Test	0.000	08:44:40
Air Blank	0.000	08:45:08
Wait	ABT*	
Air Blank	0.000	08:46:02

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
01/30/2014

08:46:29

3 of 7<sub>sc</sub>  
TSH

---

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:47:27
Diagnostic Test	Pass	08:48:02
Air Blank	0.000	08:48:29
0.100 Cal Check	0.097	08:48:50
Air Blank	0.000	08:49:19
Subject Test	RFI*	08:49:53
Air Blank	IPS**	08:50:02

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
01/30/2014

08:54:21

4 of 7 sc  
TS4

---

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:55:19
Diagnostic Test	Pass	08:55:53
Air Blank	0.000	08:56:21
0.100 Cal Check	0.097	08:56:42
Air Blank	0.000	08:57:11
Subject Test	0.174*	09:00:18
Air Blank	0.000	09:00:51
Air Blank	0.000	09:01:19
Subject Test	XXX**	09:02:05
Air Blank	0.000	09:02:38

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
01/30/2014 09:06:43

5 of 7<sub>sc</sub>  
Tbh

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:07:12
0.100 Cal Check	0.096	09:07:33
Air Blank	0.000	09:08:02
Air Blank	0.000	09:10:30

\*Input Timed Out

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 sc  
Bh  
2.12.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 1/30/14 TIME 839  
Intoxilyzer Serial # 80-001199 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
01/30/2014 08:39:22

7 of 7 sc  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:39:51
0.100 Cal Check	0.097	08:40:12
Air Blank	0.000	08:40:41

Preventative Maintenance  
Performed -- Timer Reset