



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3.11.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/26/14 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 1/30/14 and 2/26/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


3/4/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24th, day of March, 20 14.

Elma Escarcega

Notary Public

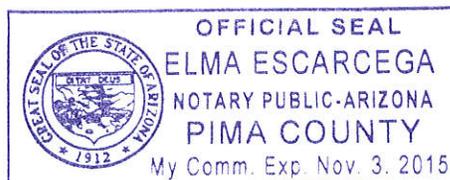

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc
TBH
3-11-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/26/14 TIME 1338

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
02/26/2014

13:40:52

2 of 7 sc
Bh

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:41:44
Diagnostic Test	Pass	13:42:19
Air Blank	0.000	13:42:47
0.100 Cal Check	0.097	13:43:07
Air Blank	0.000	13:43:37
Subject Test	0.000	13:44:03
Air Blank	0.000	13:44:30
Wait	ABT*	
Air Blank	0.000	13:45:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
02/26/2014

13:47:11

3077_{sc}
TS4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:48:07
Diagnostic Test	Pass	13:48:42
Air Blank	0.000	13:49:09
0.100 Cal Check	0.097	13:49:30
Air Blank	0.000	13:50:00
Subject Test	RFI*	13:50:25
Air Blank	IPS**	13:50:34

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
02/26/2014 13:55:54

4 of 7 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:57:03
Diagnostic Test	Pass	13:57:39
Air Blank	0.000	13:58:06
0.100 Cal Check	0.097	13:58:27
Air Blank	0.000	13:58:56
Subject Test	0.133*	14:02:03
Air Blank	0.000	14:02:35
Air Blank	0.000	14:03:02
Subject Test	XXX**	14:03:17
Air Blank	0.000	14:03:51

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
02/26/2014 14:04:22

5 of 7 sc
Bu

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:04:50
0.100 Cal Check	0.096	14:05:11
Air Blank	0.000	14:05:41

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TSN
3.11.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/26/14 TIME 1338
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

Core Version: 8105.48

02/26/2014

13:38:39

7 of 7 sc
TS4

Standard Lot#: 30412100A1

Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:39:08
0.100 Cal Check	0.096	13:39:29
Air Blank	0.000	13:39:58

Preventative Maintenance
Performed -- Timer Reset