



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY  
4.14.2014

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/27/14 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 2/26/14 and 3/27/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione  
Sarah L. Cione

4/1/14  
ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15<sup>th</sup>, day of April, 20 14.

Elma Escarcega  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

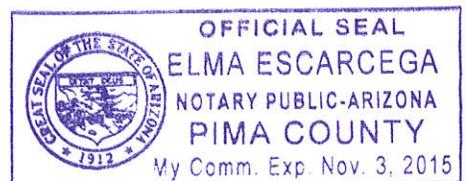


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1077<sub>sc</sub>  
TS4  
4.14.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/27/14 TIME 945

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/27/2014

09:51:48

2 of 7  
TSK

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:52:33
Diagnostic Test	Pass	09:53:08
Air Blank	0.000	09:53:35
0.100 Cal Check	0.097	09:53:56
Air Blank	0.000	09:54:25
Subject Test	0.000	09:54:51
Air Blank	0.000	09:55:18
Wait	ABT*	
Air Blank	0.000	09:56:00

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/27/2014

09:56:28

3 of 7 SC  
TBH

---

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:57:13
Diagnostic Test	Pass	09:57:48
Air Blank	0.000	09:58:15
0.100 Cal Check	RFI*	09:58:23
Air Blank	IPS**	09:58:40

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/27/2014 10:00:10

4 of 7<sub>sc</sub>  
TSh

---

Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:00:55
Diagnostic Test	Pass	10:01:31
Air Blank	0.000	10:01:58
0.100 Cal Check	0.096	10:02:19
Air Blank	0.000	10:02:48
Subject Test	0.308*	10:05:56
Air Blank	0.000	10:06:28
Air Blank	0.000	10:06:55
Subject Test	XXX**	10:07:28
Air Blank	0.000	10:08:02

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/27/2014 10:10:25

5 of 7  
x  
BY

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:10:53
0.100 Cal Check	0.096	10:11:14
Air Blank	0.000	10:11:44

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
sc  
TS4  
4.14.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 3/27/14 TIME 945  
Intoxilyzer Serial # 80-001199 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
03/27/2014 09:45:44

7 of 7  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. SULLIVAN #34971  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:46:12
0.100 Cal Check	0.096	09:46:33
Air Blank	0.000	09:47:02

Preventative Maintenance  
Performed -- Timer Reset