



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 9.29.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/11/14 Cal 1 0.095 AC Cal 2 0.095 AC

Between the dates of 8/4/14 and 9/11/14 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione
 Sarah L. Cione

9/15/14
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15th, day of September, 20 14.

Elma Escarcega
 Notary Public

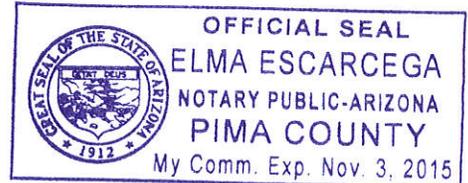

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
5c
154
9.29.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/11/14 TIME 1007

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
09/11/2014

10:09:25

2 of 8 sc

PH

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:10:20
Diagnostic Test	Pass	10:10:55
Air Blank	0.000	10:11:22
0.100 Cal Check	0.094	10:11:43
Air Blank	0.000	10:12:13
Subject Test	0.000	10:12:52
Air Blank	0.000	10:13:20
Wait	ABT*	
Air Blank	0.000	10:14:33

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
09/11/2014 10:16:46

3 of 8
SL
TJ

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:17:41
Diagnostic Test	Pass	10:18:16
Air Blank	0.000	10:18:43
0.100 Cal Check	0.095	10:19:04
Air Blank	0.000	10:19:33
Subject Test	RFI*	10:20:02
Air Blank	IPS**	10:20:13

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
09/11/2014 10:21:03

4 of 8 SC
PS4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:21:51
Diagnostic Test	Pass	10:22:26
Air Blank	0.000	10:22:54
0.100 Cal Check	0.095	10:23:14
Air Blank	0.000	10:23:44
Subject Test	XXX*	10:24:34
Air Blank	0.000	10:25:06

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
09/11/2014 10:25:12

5078_{SC}
TB4

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:25:52
Diagnostic Test	Pass	10:26:27
Air Blank	0.000	10:26:55
0.100 Cal Check	0.095	10:27:16
Air Blank	0.000	10:27:45
Subject Test	0.109*	10:30:53
Air Blank	0.000	10:31:24
Air Blank	0.000	10:31:52
Subject Test	ABT**	10:32:42
Air Blank	0.000	10:33:10

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
09/11/2014 10:34:18

6 of 8
SC
104

Standard Lot#: 30412100A1
Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:34:46
0.100 Cal Check	0.095	10:35:07
Air Blank	0.000	10:35:37

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8 SC
134
9.29.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/11/14 TIME 1007
Intoxilyzer Serial # 80-001199 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

Core Version: 8105.48

09/11/2014

10:07:00

Standard Lot#: 30412100A1

Last Changed By: S. SULLIVAN #34971

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test g/210L Time

Air Blank 0.000 10:07:29
0.100 Cal Check 0.095 10:07:49
Air Blank 0.000 10:08:19

Preventative Maintenance
Performed -- Timer Reset

80/8_{sc}

104