



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

Handwritten initials and date: 2/3/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-26-2015 Cal 1 0.100 AC Cal 2 — AC

Between the dates of — and — the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 3-26-2015.
The instrument was returned to service on —.

Handwritten signature of Terry S. Gallegos

Terry S. Gallegos


1-27-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 27th, day of January, 2016.

Handwritten signature of Elma Escarcega

Notary Public

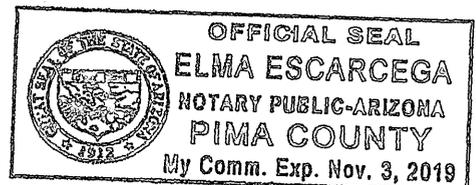

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

1077
TS4
JL
2/3/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3-26-2015 TIME 1220
Intoxilyzer Serial # 80-001199 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START" -
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed.
- () 5. Calibration check completed. Test results 0. _____ AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

EEProm failure Standby mode

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Could not complete the check list. I took the instrument to the Crime Lab for further testing. TS4

At 17:16 I tried a subject test: EEPROM checksum test Fail TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/26/2015 17:16:17

2 of 7
734
All

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:17:44
Diagnostic Test	Fail*	17:17:49
Air Blank	0.000	17:18:17

* EEPROM Checksum Test

Not a Successfully Completed Test Sequence

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

3 of 7
TSY
JDA

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3-26-2015 TIME 1719
Intoxilyzer Serial # 80-001199 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER. *Used Technician card to go to control test menu.*
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/26/2015 17:19:11

4 of 7
TS4
DB

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:19:40
0.100 Cal Check	0.100	17:20:01
Air Blank	0.000	17:20:31

Preventative Maintenance
Performed -- Timer Reset

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
Core Version: 8105.48
03/26/2015 17:21:07

5 of 7
TS4
AD

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:21:56
Diagnostic Test	Fail*	17:22:01
Air Blank	0.000	17:22:29

* EEPROM Checksum Test

Tried another subject test with Diagnostic Fail.
TS4

Not a Successfully Completed Test Sequence



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

60 of 7
TBH
[Signature]

Date: 3.26.2015

QAS: Terry Gallegos

Intoxilyzer Serial Number: 80-001199

Reason For Service: Diagnostic Fail during subject test.
EEPROM checksum test message. I downloaded all records
to COBRA. I tried the Flash update from COBRA software
and EEPROM checksum test message fail.

Date Out of Service: _____

Calibration Check completed prior to removal from service

Date: 3.26.2015

SQAP completed prior to removal from service

Date: No

Corrective Action: Sent to CMI, Inc.

Sent to CMI on

The instrument was stored at TPD Crime Lab since
March 26, 2015. TBH



Service Evaluation Form

7077
Tg
JL

This form MUST be completed and enclosed with instrument to be serviced.

Note: Please ship items in their original shipping container.

Contact information: Customer Number 105075 (contact Customer Service)
Name TERRY GALLEGOS Phone: (520) 837-7677
Fax: (520) 791-4642 Email: TERRY.GALLEGOS@TUCSONAZ.GOV

Bill to Address:
FINANCE/ACCOUNTS PAYABLE
PO BOX 27450
TUCSON, AZ 85726-7450

Ship to Address:
TUCSON POLICE DEPARTMENT CRIMELAB
1306 W. MIRAGE MILE
TUCSON, AZ 85705

Instrument Serial Number: 80-001199

Detailed Description of Problem:
Diagnostic fail prior to subject test. EEPROM checksum test message.
Flash update from COBRA software did not clear EEPROM checksum
test message fail.

***Hazardous Material Warning! - DO NOT return gas cylinders with instrument! ***

Note: An evaluation fee of \$39.50 for handheld instruments or \$79.00 for IR evidential instruments will apply to estimates which are not repaired.

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
Purchase Order Number 37768

Authorized By:

TERRY GALLEGOS TOXICOLOGY COORDINATOR Ship to:

Name (Please Print)

Title

Terry Gallegos
Signature

1-26-2016
Date

CMI, Inc.
Attn: Service Dept.
316 East Ninth Street
Owensboro, KY 42303

No, please send estimate before repairs are made. An estimate will be faxed before performing any repairs and may cause delays in service.