



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL  
6/29/16

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-16-2016 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 5-26-2016 and 6-16-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Terry S. Gallegos

6-29-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 29th, day of June, 2016.

Notary Public

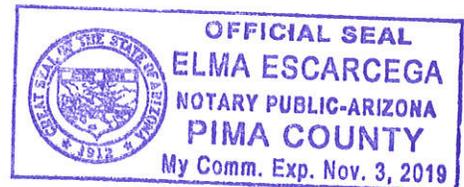

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1087  
TS4  
WLL  
6/29/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6.16.2016 TIME 1324  
INTOXILYZER SERIAL # 80-001199 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
06/16/2016 13:29:07

2 of 7  
By  
WLL

Standard Lot#: 15615100A1  
Last Changed By: W. LOOMIS #53545

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:29:53
Diagnostic Test	Pass	13:30:28
Air Blank	0.000	13:30:56
0.100 Cal Check	0.098	13:31:16
Air Blank	0.000	13:31:46
Subject Test	0.000	13:32:14
Air Blank	0.000	13:32:41
Wait	ABT*	
Air Blank	0.000	13:33:15

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
06/16/2016

13:51:35

3087  
TS4  
WLL

---

Standard Lot#: 15615100A1  
Last Changed By: W. LOOMIS #53545

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:52:35
Diagnostic Test	Pass	13:53:10
Air Blank	0.000	13:53:38
0.100 Cal Check	RFI*	13:53:54
Air Blank	IPS**	13:54:08

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
06/16/2016

13:55:40

4 of 7  
TS  
WLL

---

Standard Lot#: 15615100A1  
Last Changed By: W. LOOMIS #53545

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:56:27
Diagnostic Test	Pass	13:57:02
Air Blank	0.000	13:57:29
0.100 Cal Check	0.098	13:57:50
Air Blank	0.000	13:58:19
Subject Test	0.095*	14:01:27
Air Blank	0.000	14:01:58
Air Blank	0.000	14:02:26
Subject Test	XXX**	14:03:20
Air Blank	0.000	14:03:54

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
06/16/2016 14:04:16

Soft  
By  
WLL

-----  
Standard Lot#: 15615100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:04:45
0.100 Cal Check	0.097	14:05:05
Air Blank	0.000	14:05:35

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
B4  
WLL  
6/29/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6-16-2016 TIME 1326  
Intoxilyzer Serial # 80-001199 Location DUI-HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001199  
Core Version: 8105.48  
06/16/2016

13:26:50

707  
154  
WLL

-----  
Standard Lot#: 15615100A1  
Last Changed By: W. LOOMIS #53545  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:27:32
0.100 Cal Check	0.099	13:27:53
Air Blank	0.000	13:28:23

Preventative Maintenance  
Performed -- Timer Reset