



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

je

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-28-2006 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 12-1-2006 and 12-28-2006 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry Gallegos
Terry S. Gallegos

1.4.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

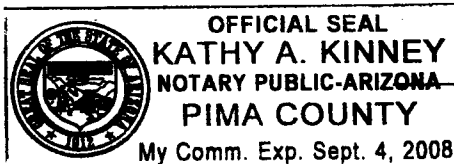

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police DepartmentDATE 12-28-2006 TIME 1614INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed ✓✓
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 097
- AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

2 of 8
TAL

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
12/28/2006

16:16:20

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:17:15
Diagnostic Test	Pass	16:17:42
Air Blank	0.000	16:18:10
0.100 Cal Check	0.098	16:18:28
Air Blank	0.000	16:18:58
Subject Test	0.000	16:19:54
Air Blank	0.000	16:20:21
Five Minute Wait	ABT*	
Air Blank	0.000	16:21:07

*Sequence Aborted

Not a Successfully Completed Test Sequence

3048
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
12/28/2006

ju

16:22:01

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:22:53
Diagnostic Test	Pass	16:23:20
Air Blank	0.000	16:23:47
0.100 Cal Check	0.098	16:24:06
Air Blank	0.000	16:24:36
Subject Test	RFI*	16:24:59
Air Blank	0.000	16:25:27

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002201

12/28/2006

16:25:34

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	16:26:10
Air Blank	0.000	16:26:38

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
12/28/2006

16:27:04

50f8
784

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:27:55
Diagnostic Test	Pass	16:28:23
Air Blank	0.000	16:28:50
0.100 Cal Check	0.099	16:29:09
Air Blank	0.000	16:29:38
Subject Test	0.015*	16:32:47
Air Blank	0.000	16:33:16
Air Blank	0.000	16:33:44
Subject Test	XXX**	16:34:10
Air Blank	0.000	16:34:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

12/28/2006

16:35:04

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:35:31
Cal Check	0.097	16:35:50
Air Blank	0.000	16:36:19


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8
184

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

ye

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28047 AGENCY Tucson Police Department
DATE 12-28-2006 TIME 1614
Intoxilyzer Serial # 80-002201 Location DVI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

12/28/2006

16:14:22

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:14:49
Cal Check	0.098	16:15:08
Air Blank	0.000	16:15:38


Operator's Signature