



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

ms

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-26-2007 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 12-28-2006 and 1-26-2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos
Terry S. Gallegos

2-27-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 27, day of

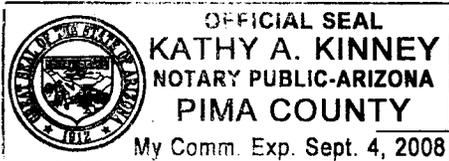

February, 2007.

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

108 11 7ms
TBY

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-26-2007 TIME 1704
INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
01/26/2007 17:06:17

2 of 11
BY

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:07:13
Diagnostic Test	Pass	17:07:41
Air Blank	0.000	17:08:08
0.100 Cal Check	0.098	17:08:27
Air Blank	0.000	17:08:56
Subject Test	0.000	17:09:31
Air Blank	0.000	17:09:58
Five Minute Wait	ABT*	
Air Blank	IPS**	17:10:26

*Sequence Aborted
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
01/26/2007 17:10:39

3 of 11
176

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:11:43
Diagnostic Test	Pass	17:12:10
Air Blank	0.000	17:12:38
0.100 Cal Check	INT*	17:12:58
Air Blank	AMB**	17:13:29

*Interferent Detect
**Ambient Fail

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
01/26/2007 17:14:03

40f 11
TSY

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:15:02
Diagnostic Test	Pass	17:15:30
Air Blank	0.000	17:15:57
0.100 Cal Check	0.097	17:16:16
Air Blank	RFI*	17:16:33
Air Blank	0.000	17:17:01

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
01/26/2007 17:17:37

508 11
124

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:18:34
Diagnostic Test	Pass	17:19:01
Air Blank	0.000	17:19:29
0.100 Cal Check	0.097	17:19:48
Air Blank	0.000	17:20:17
Subject Test	0.413*	17:23:25
Air Blank	0.000	17:24:00
Air Blank	0.000	17:24:27
Subject Test	XXX**	17:24:48
Air Blank	0.000	17:25:22

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
01/26/2007 17:25:58

60f 11
124 23

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:26:25
Cal Check	0.097	17:26:44
Air Blank	0.000	17:27:14


Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

*70f 11
104*

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 2884 AGENCY Tucson Police Department
DATE 1-26-2007 TIME 1704
Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000 .
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
01/26/2007 17:04:15

808 11
TSH

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:04:42
Cal Check	0.097	17:05:01
Air Blank	0.000	17:05:31


Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

907 11
TSC

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 1.26.2007
 QAS: Terry Gallegos 28847
 Intoxilyzer Serial Number: 80-002201
 New Lot Number: 627101G
 New Expiration Date: 10.03.2008

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000 .
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
01/26/2007 17:29:23

100f 11
TM

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:29:50
Cal Check	0.099	17:30:09
Air Blank	0.000	17:30:38



Operator's Signature

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002201
01/26/2007 17:37:34

11 of 11
TJ4

Standard Lot#: 627101G
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:38:29
Diagnostic Test	Pass	17:38:56
Air Blank	0.000	17:39:24
0.100 Cal Check	0.099	17:39:43
Air Blank	0.000	17:40:12
Subject Test	RFI*	17:40:30
Air Blank	0.000	17:40:58

*RFI Detect

Tried to get RFI on the Cal Check,
but it did not register during that part
of the sequence. BJ

Not a Successfully Completed Test Sequence