



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

2ms  
7-13-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-21-2007 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 5-24-2007 and 6-21-2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.


Terry S. Gallegos

Terry S. Gallegos

7-13-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13, day of July, 2007.

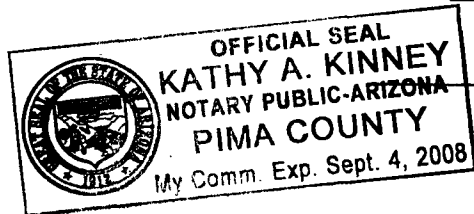

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 80001 of 8  
TS4 #  
71307

## STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police DepartmentDATE 6-21-2007 TIME 1741INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- 
- (✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
- 
- (✓) 2. Error recognition logic system functioning.
- 
- Not a Successfully Completed Test Sequence printed ✓✓✓
- 
- (✓) 3. Proper sample recognition system.
- 
- Not a Successfully Completed Test Sequence printed
- 
- Deficient sample printed.
- 
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
- 
- Result: 0.
- 099
- AC

Instrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

---

---

---

---

SIGNATURE



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002201

06/21/2007

17:42:51

20f8  
TS4 <sup>mm</sup> 7/3/07

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:43:46
Diagnostic Test	Pass	17:44:14
Air Blank	0.000	17:44:41
0.100 Cal Check	0.100	17:45:00
Air Blank	0.000	17:45:29
Subject Test	0.000	17:46:10
Air Blank	0.000	17:46:38
Five Minute Wait	ABT*	
Air Blank	0.000	17:47:29

\*Sequence Aborted

### **Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002201  
06/21/2007

17:48:20

30f8  
TSy MB  
7/30/07

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:49:08
Diagnostic Test	Pass	17:49:36
Air Blank	0.000	17:50:03
0.100 Cal Check	INT*	17:50:22
Air Blank	RFI**	17:50:30

\*Interferent Detect

\*\*RFI Detect

*TSy reference, just before the*  
*I keyed the radio during the Cal Check, and it*  
*caused the interferent detect.*

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002201  
06/21/2007

17:51:00

40F8  
TS4 7/13/07

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:51:49
Diagnostic Test	Pass	17:52:16
Air Blank	0.000	17:52:44
0.100 Cal Check	0.101	17:53:02
Air Blank	0.000	17:53:32
Subject Test	0.297*	17:56:41
Air Blank	0.000	17:57:14
Air Blank	0.000	17:57:41
Subject Test	XXX**	17:58:09
Air Blank	IPS***	17:58:18

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

\*\*\*Improper Sample

### **Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

06/21/2007

17:58:56

50f8  
TSy m  
71507

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	AMB*	17:59:24
Air Blank	0.000	17:59:53

\*Ambient Fail

  
Operator's Signature

I was close to the instrument during the air blank  
and just after the Invalid Sample test; I  
caused the ambient fail.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

06/21/2007

17:59:58

60f8  
B4 mm  
7/2/07

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:00:25
Cal Check	0.099	18:00:44
Air Blank	0.000	18:01:13

  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6.21.2007 TIME 1741  
Intoxilyzer Serial # 80-002201 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002201

06/21/2007

17:41:07

80f8

TS4 MS  
7/3/07

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	17:41:34
Cal Check	0.100	17:41:53
Air Blank	0.000	17:42:23

  
Operator's Signature