



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

*MMS
10-25-07*

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.19.2007 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8.17.2007 and 9.19.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos
Terry S. Gallegos

10.25.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25, day of October, 2007.

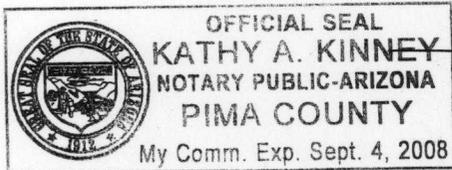

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

MMR
10-25-07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 8

STANDARD QUALITY ASSURANCE PROCEDURE

TJH

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9-19-2007 TIME 1734

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓✓
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Ex G-2 performed at 33 days. TJH

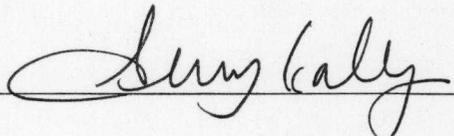SIGNATURE 

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:35:54

20f8
TS4
7ms
10-25-07

Standard Lot#: 706602I
Last Changed By: J. OLIUARES #36552

OAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:36:42
Diagnostic Test	Pass	17:37:11
Air Blank	0.000	17:37:38
0.100 Cal Check	0.099	17:37:57
Air Blank	0.000	17:38:27
Subject Test	0.000	17:38:43
Air Blank	0.000	17:39:10
Five Minute Wait	ABT*	
Air Blank	0.000	17:39:44

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

30f8
TS4
mb
10-2507

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:39:55

Standard Lot#: 7066021
Last Changed By: J. OLIUARES #36552

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:40:46
Diagnostic Test	Pass	17:41:15
Air Blank	0.000	17:41:42
0.100 Cal Check	0.099	17:42:01
Air Blank	RFI*	17:42:24
Air Blank	IPS**	17:42:37

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

40f8
TJ
MB
10-25-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:42:51

Standard Lot#: 7066021
Last Changed By: J. OLIUARES #36552

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	AMB*	17:43:43
Air Blank	0.000	17:44:12

*Ambient Fail

Not a Successfully
Completed
Test Sequence

5048
TS4
2/15
10-25-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:44:29

Standard Lot#: 7066021
Last Changed By: J. OLIUARES #36552

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQA, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:45:27
Diagnostic Test	Pass	17:45:56
Air Blank	0.000	17:46:23
0.100 Cal Check	0.099	17:46:42
Air Blank	0.000	17:47:12
Subject Test	0.160*	17:50:20
Air Blank	0.000	17:50:53
Air Blank	0.000	17:51:20
Subject Test	XXX**	17:51:46
Air Blank	0.000	17:52:20

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence


608
TS4
7/16
10-25-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:52:57

Standard Lot#: 7066021
Last Changed By: J. OLIVARES #36552

DAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:53:25
Cal Check	0.099	17:53:43
Air Blank	0.000	17:54:13

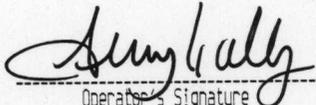

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7078
TS4
WMS
10-25-07

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9-19-2007 TIME 1734
Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

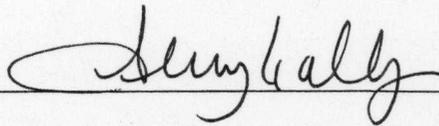SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

80f8
TJh
10-25-07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002201
09/19/2007 17:34:10

Standard Lot#: 7066021
Last Changed By: J. OLIVARES #36552

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:34:38
Cal Check	0.099	17:34:57
Air Blank	0.000	17:35:26



Operator's Signature
1