



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

CU  
2-8-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: - Cal 1 N/A AC Cal 2 N/A AC

Between the dates of - and - the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 11/2/2007.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/8/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 8, day of February, 2008.

Kathy A Kinney

Notary Public

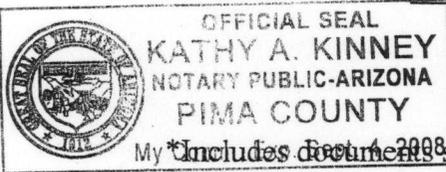

My Commission Expires 2/28/2008 as specified in A.A.C. R13-10-104.A.5



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

CU  
2-8-08

INTOXILYZER SERVICE RECORD

1051  
B4

Date: 11-8-2007

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-002201

Reason For Service: Fails "modem test" during diagnostics.

Date Out of Service: 11-02-2007

Calibration Check completed prior to removal from service

Date: NO

SQAP completed prior to removal from service

Date: NO

Corrective Action:  Sent to CMI, Inc.

COBRA Data download prior to shipping, 11-8-2007



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

CM  
2-8-08

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

1 of 2  
PJ

Date: 11-02-07  
QAS: OLIVARES, JOSE #36552  
Intoxilyzer Serial Number: 80-002201  
New Lot Number: 719802 I  
New Expiration Date: 23 JUL 2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
  - Enter Target Value: 0.100
  - Enter Control Lot #: From the tank label
  - Enter Exp Date: From the tank label (mm/dd/yyyy)
  - Edit info Y/N? Y or N
  - Update Std Info? Y ENTER
  - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Jose Olivares #36552

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

# Tucson Police Department Crime Laboratory


CM  
2-8-08  
2 of 2  
154

INTOXLYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002201  
11/02/2007 01:59:00

Standard Lot#: 7198021  
Last Changed By: J. OLIVARES #36552

QAS: J. OLIVARES #36552  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	01:59:28
Cal Check	0.100	01:59:47
Air Blank	0.000	02:00:16

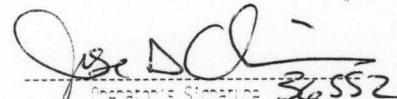  
Operator's Signature 36552