



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TS4  
5-12-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/9/08 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 4/10/08 and 5/9/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione


Sarah L. Cione

5/12/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 12, day of May, 2008.

Kathy A Kinney  
Notary Public

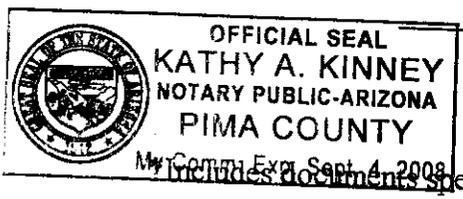

My Comm. Exp. Sept. 4, 2008  
These documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 9 sc

TS4  
5.12.08

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/9/08 TIME 9:51

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
05/09/2008

09:53:26

2 of 9 sc

BH

Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:54:15
Diagnostic Test	Pass	09:54:42
Air Blank	0.000	09:55:09
0.100 Cal Check	0.099	09:55:28
Air Blank	0.000	09:55:58
Subject Test	0.000	09:56:36
Air Blank	0.000	09:57:04
Five Minute Wait	ABT*	
Air Blank	0.000	09:57:51

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
05/09/2008

09:57:56

3 of 9 sc  
SM

Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:58:46
Diagnostic Test	Pass	09:59:14
Air Blank	RFI*	09:59:24
Air Blank	0.000	09:59:52

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
05/09/2008

10:00:40

4 of 9 sc  
BH

Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	10:01:55
Air Blank	0.000	10:02:24

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
05/09/2008

10:02:28

5 of 9 sc

B4

Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:03:52
Diagnostic Test	Pass	10:04:20
Air Blank	0.000	10:04:47
0.100 Cal Check	0.099	10:05:06
Air Blank	0.000	10:05:36
Subject Test	0.067*	10:08:44
Air Blank	0.000	10:09:16
Air Blank	0.000	10:09:44
Subject Test	RGE**	10:10:09
Air Blank	0.000	10:10:45

\*Deficient Sample

\*\*Range Exceeded

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
05/09/2008

10:12:16

6 of 9 sc

By

Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:13:01
Diagnostic Test	Pass	10:13:28
Air Blank	0.000	10:13:56
0.100 Cal Check	0.099	10:14:14
Air Blank	0.000	10:14:44
Subject Test	XXX*	10:15:02
Air Blank	0.000	10:15:36

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
05/09/2008

10:15:54

7 of 9 SC

-----  
Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552  
-----

164

-----  
QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:16:22
Cal Check	0.099	10:16:41
Air Blank	0.000	10:17:10

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

8 of 9 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
5.12.08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/9/08 TIME 9:51  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
05/09/2008

09:51:33

-----  
Standard Lot#: 732503I  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:52:00
Cal Check	0.099	09:52:19
Air Blank	0.000	09:52:49

*Sarah Cione*

Operator's Signature

9 of 95c

104