



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-10-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/5/08 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 5/9/08 and 6/5/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/9/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9, day of June, 2008.

Kathy A Kinney

Notary Public

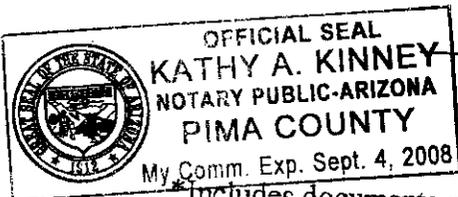

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

B4
6-10-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/5/08 TIME 9:15

INTOXILYZER SERIAL # 80-002201 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

134

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/05/2008 09:17:50

Standard Lot#: 7325031
Last Changed By: J. OLIVARES #36552

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:18:47
Diagnostic Test	Pass	09:19:15
Air Blank	0.000	09:19:43
0.100 Cal Check	0.098	09:20:02
Air Blank	0.000	09:20:31
Subject Test	0.000	09:20:48
Air Blank	0.000	09:21:15
Five Minute Wait	ABT*	
Air Blank	0.000	09:21:51

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/05/2006 09:22:10

Standard Lot#: 7325031
Last Changed By: J. OLIVARES #36552

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:23:07
Diagnostic Test	Pass	09:23:36
Air Blank	RFI*	09:23:48
Air Blank	IPS**	09:23:57

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7 sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/05/2008 09:24:15

Standard Lot#: 7325031
Last Changed By: J. OLIVARES #36552

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:25:12
Diagnostic Test	Pass	09:25:41
Air Blank	0.000	09:26:08
0.100 Cal Check	0.097	09:26:27
Air Blank	0.000	09:26:57
Subject Test	0.022*	09:30:05
Air Blank	0.000	09:30:37
Air Blank	0.000	09:31:05
Subject Test	XXX**	09:31:29
Air Blank	0.000	09:32:03

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7 sc

ph

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002201
06/05/2008 09:35:35

Standard Lot#: 732503T
Last Changed By: J. OLIVARES #36552

OAS: S. CIONE #50616
TUCSON PD

Test g/210L Time

Air Blank 0.000 09:36:02
Cal Check 0.097 09:36:21
Air Blank 0.000 09:36:50

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

By
6-10-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/5/08 TIME 9:15
Intoxilyzer Serial # 80-002201 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 75c
BH

INTOXILYZER 8600
Location: ARIZONA
Serial Number: 80-002201
06/05/2008 09:15:49

Standard Lot#: 7325031
Last Changed By: J. OLIVARES #36552

GAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:16:16
Cal Check	0.098	09:16:35
Air Blank	0.000	09:17:05

Sarah Cione
Operator's Signature