



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

*u*  
7-10-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002201 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/3/08 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 6/5/08 and 7/3/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

7/3/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3, day of July, 2008.

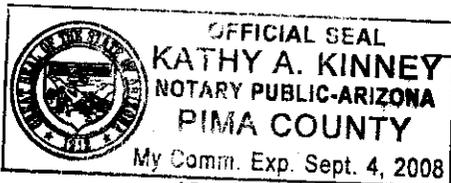

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 sc  
CM  
7-10-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/3/08 TIME 12:43

INTOXILYZER SERIAL # 80-002201 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7 sc  
cu

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/03/2008 12:44:55

---

Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:45:41
Diagnostic Test	Pass	12:46:09
Air Blank	0.000	12:46:36
0.100 Cal Check	0.099	12:46:55
Air Blank	0.000	12:47:25
Subject Test	0.000	12:47:45
Air Blank	0.000	12:48:12
Five Minute Wait	ABT*	
Air Blank	0.000	12:48:46

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/03/2008

12:48:51

3 of 7 sc  
lm

---

Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:49:36
Diagnostic Test	Pass	12:50:04
Air Blank	RFI*	12:50:14
Air Blank	IPS**	12:50:23

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 sc

Ch

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002201  
07/03/2008 12:50:33

Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:51:14
Diagnostic Test	Pass	12:51:42
Air Blank	0.000	12:52:09
0.100 Cal Check	0.098	12:52:28
Air Blank	0.000	12:52:58
Subject Test	0.286*	12:56:06
Air Blank	0.000	12:56:39
Air Blank	0.000	12:57:07
Subject Test	XXX**	12:57:34
Air Blank	0.000	12:58:08

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

OK

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
07/03/2008

12:58:36

-----  
Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:59:03
Cal Check	0.098	12:59:22
Air Blank	0.000	12:59:52

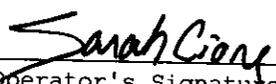
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

CU

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

7-10-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 7/3/08 TIME 12:43  
Intoxilyzer Serial # 80-002201 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002201  
07/03/2008

12:43:34

7 of 75c  
Ch

-----  
Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:44:01
Cal Check	0.099	12:44:20
Air Blank	0.000	12:44:50

Sarah Cione  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2 sc


CU

7-10-08

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 062908  
 QAS: OLIVARES, JOSE A.  
 Intoxilyzer Serial Number: 80-002201  
 New Lot Number: 735401G  
 New Expiration Date: 12-28-09

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

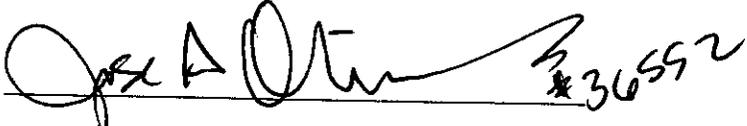SIGNATURE  #36552

2 of 2 sc

Oh

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 88-002201  
06/29/2008 01:33:09

-----  
Standard Lot#: 735401G  
Last Changed By: J. OLIVARES #36552

-----  
OAS: J. OLIVARES #36552  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	01:33:37
Cal Check	0.098	01:33:55
Air Blank	0.000	01:34:25

  
Operator's Signature 36552